
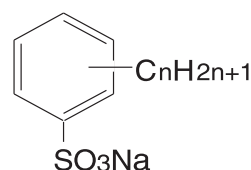


陰イオン界面活性剤
Anionic surfactants

N040324D

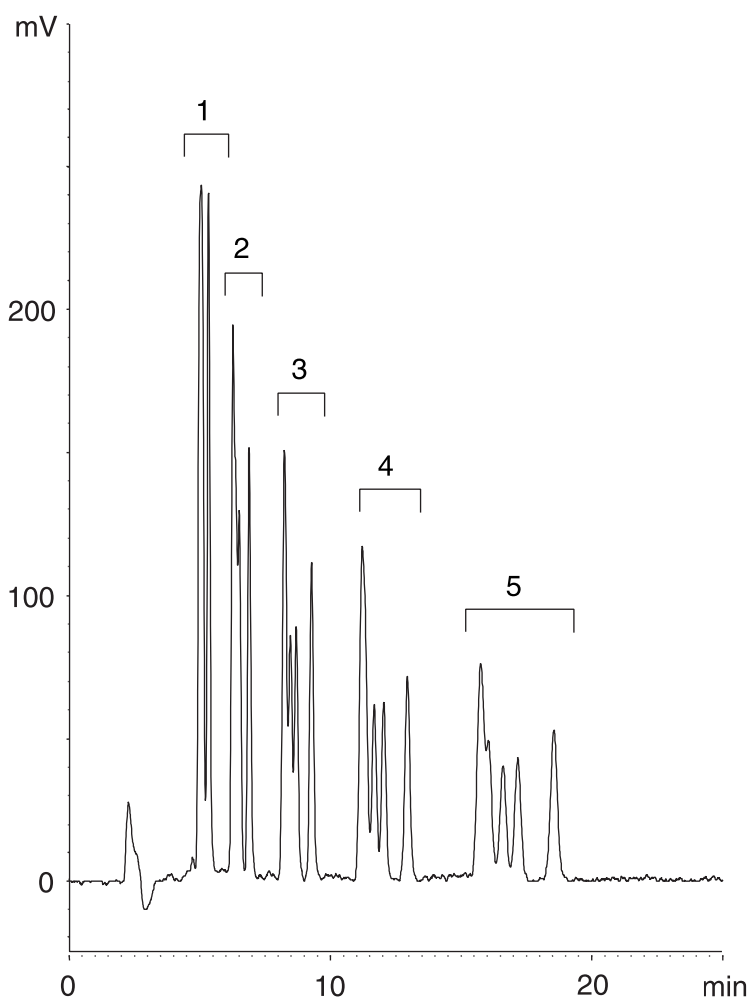

1. Sodium decylbenzenesulfonate (n=10)
2. Sodium undecylbenzenesulfonate (n=11)
3. Sodium dodecylbenzenesulfonate (n=12)
4. Sodium tridecylbenzenesulfonate (n=13)
5. Sodium tetradecylbenzenesulfonate (n=14)

Column : YMC-Pack Pro C18 (5 μ m, 12nm)
250 \times 4.6 mm I.D.
Eluent : acetonitrile/water (65/35) containing 0.1M NaClO₄
Flow rate : 1.0 mL/min
Temperature : 40 °C
Detection : FLS at Ex 221 nm, Em 284 nm
Injection : 20 μ L (2.0 μ g/mL)
平成15年改正水質基準の検査方法準拠