


殺菌薬
Disinfectants

S930617D

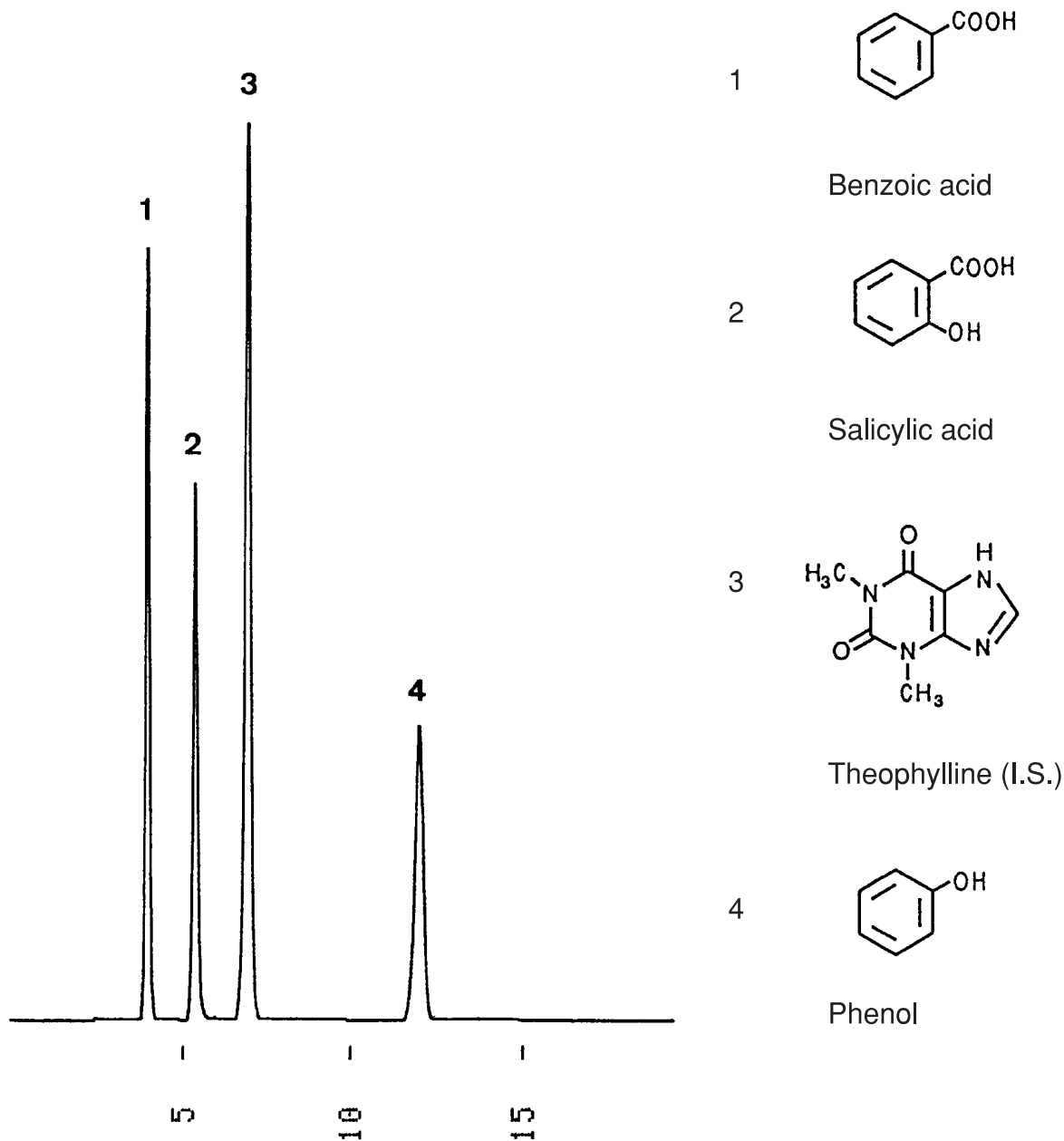

Column : YMC-Pack ODS-AL (5 μ m, 12nm)
150 \times 4.6 mm I.D.

Eluent : 100 mM KH₂PO₄-Na₂HPO₄ (pH 7.0)/methanol(75/25)

Flow rate : 0.8 mL/min

Temperature : 30 $^{\circ}$ C

Detection : UV at 270 nm, 0.13 AUFS

Injection : 3 μ L (0.1~0.8 mg/mL)

日本薬局方第12改正記載条件