



アシクロビル(シロップ・注射液)
Aciclovir syrup and injection

R100715Q

Standard solution for aciclovir syrup

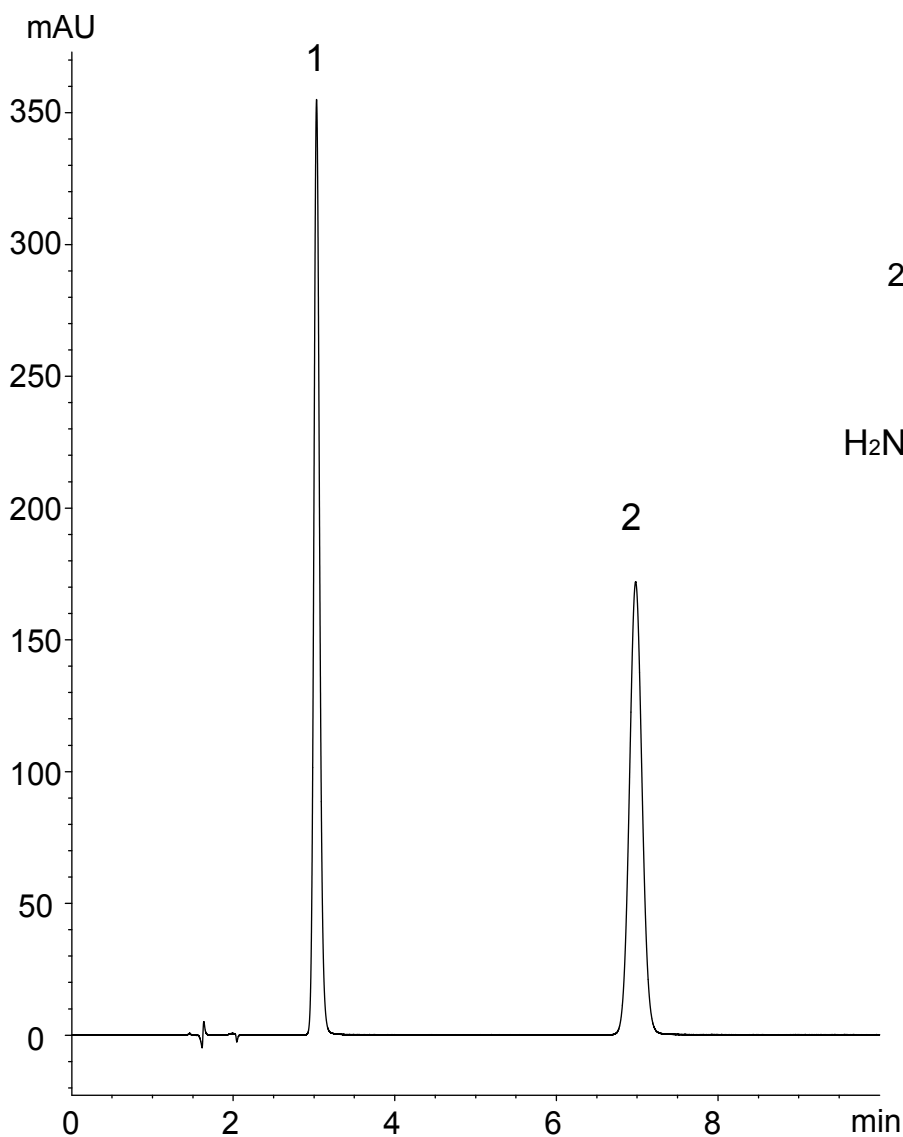


1.

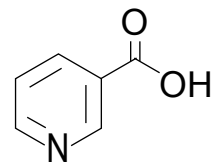

Nicotinic acid

2.

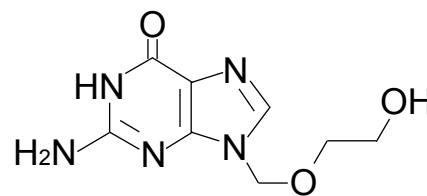

Aciclovir

Column : YMC-Triart C18 (5 μ m, 12 nm)
150 X 4.6 mm I.D.
Eluent : phosphate buffer*/methanol (950/50)
*Dissolve 1.45 g of H₃PO₄ and 25 mL of 1 mol/L CH₃COOH in water to make 900 mL → adjust pH 2.5 by 1 mol/L NaOH → add water to make 1000 mL
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 254 nm
Injection : 20 μ L (0.05 mg/mL, 0.032 mg/mL)

(日本薬局方収載原案記載条件)