
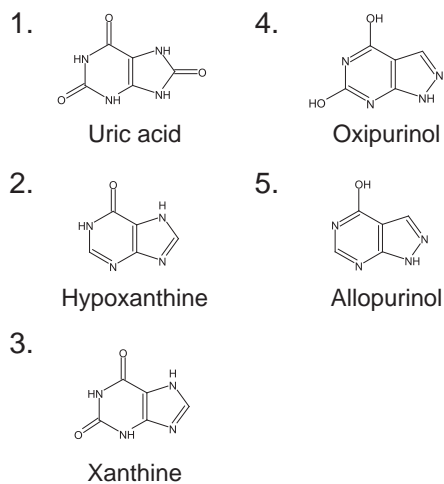


血清中のキサンチンオキシダーゼ阻害剤代謝物


Xanthine oxidase inhibitor and its metabolites in human serum

F070906A

A) Standard (100 mg/mL)

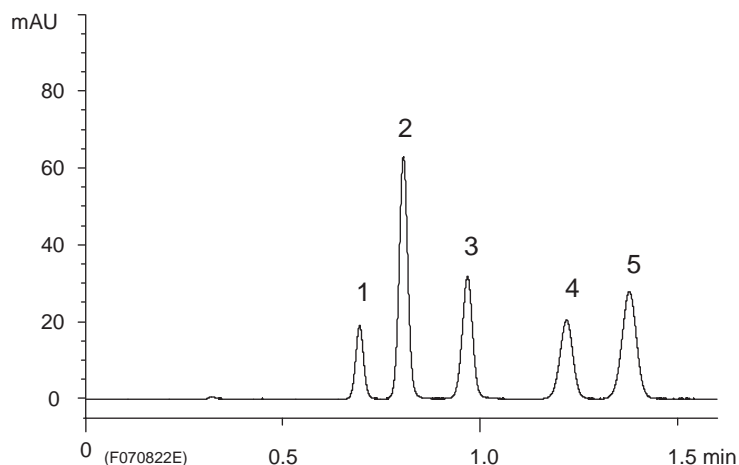


B) Spiked human serum (10 mg/mL)

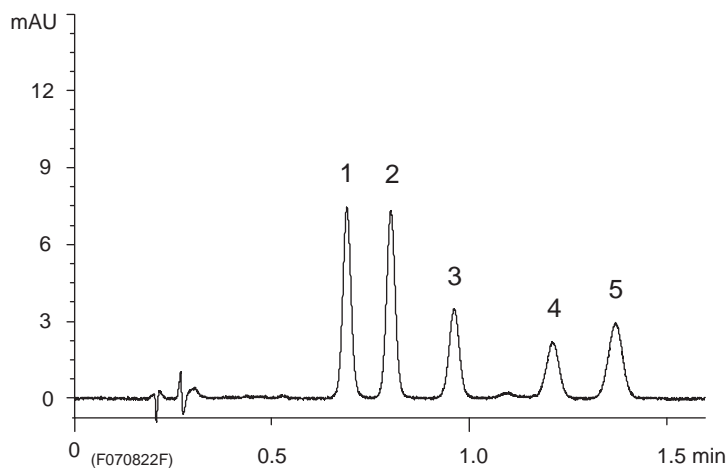

< Solid-phase extraction method >

YMC Dispo SPE C18 100 mg/1 mL

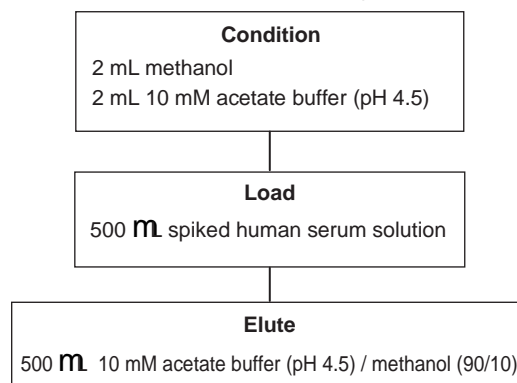

Column : YMC-UltraHT Hydrosphere C18 (2 mm, 12 nm)
50 X 2.0 mm I.D.
Eluent : 10 mM CH₃COOH-CH₃COONH₄ (pH 4.5)
Flow rate : 0.6 mL/min
Temperature : 37 °C
Detection : UV at 254 nm
Injection : 0.5 μ L