


抗生物質
Antibiotics

G910523M

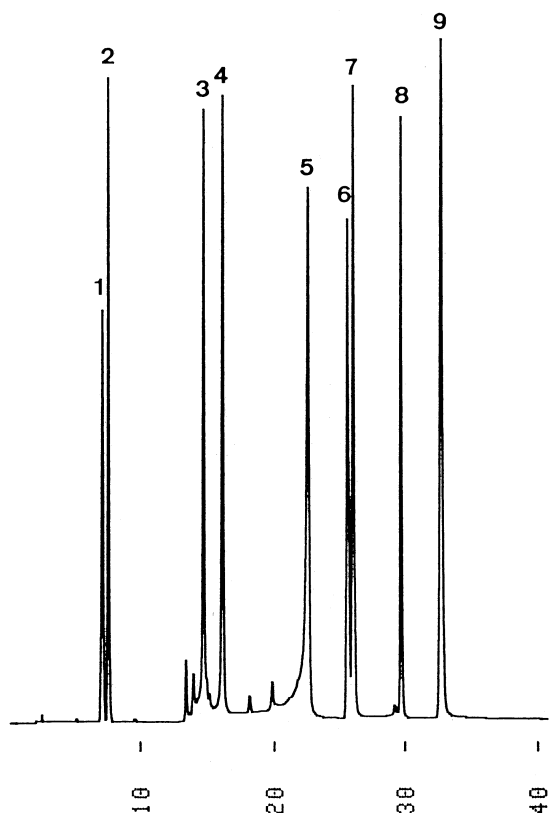

- 1.Ampicilline
- 2.Cephalexin
- 3.Oxytetracycline hydrochloride
- 4.Tetracycline
- 5.Chlorotetracycline
- 6.Chloramphenicol
- 7.Oxolinic acid
- 8.Doxorubicin hydrochloride
- 9.Nalidixic acid

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.

Eluent : A)acetonitrile/50mM NH₄H₂PO₄ (10/90)
B)acetonitrile/50mM NH₄H₂PO₄ (50/50)
0-55%B (0-30min,linear), 55%B (30-40min)

Flow rate : 1.0 mL/min

Temperature : 30 $^{\circ}$ C

Detection : UV at 250 nm, 0.16 AUFS

Injection : 15 μ L (0.01~0.15 mg/mL)