
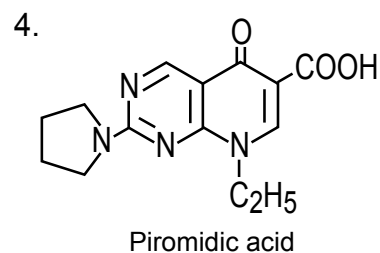
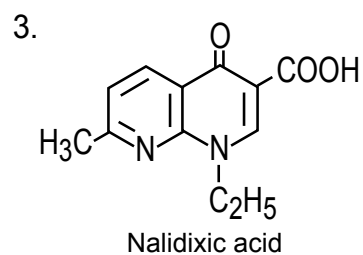
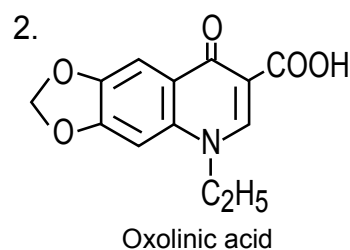
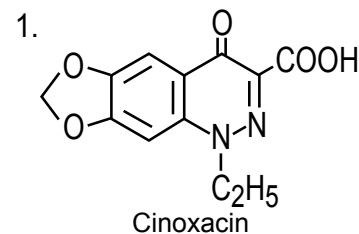


ナリジクス酸系抗生物質  
Nalidixic acid antibiotics

A010324A

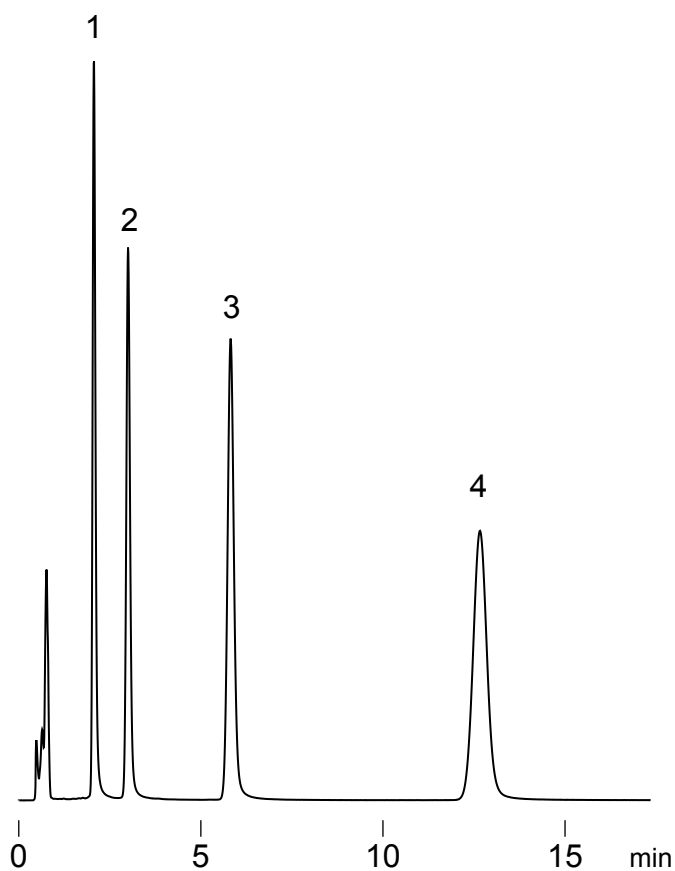

Column	: HydrosphereC18 (3 $\mu$ m, 12 nm) 50 X 2.0 mm I.D.
Eluent	: 10 mM CH <sub>3</sub> COOH-CH <sub>3</sub> COONH <sub>4</sub> (pH 4.3)/acetonitrile (75/25)
Flow rate	: 0.2 mL/min
Temperature	: 37°C
Detection	: UV at 260 nm, 0.128 AUFS
Injection	: 1 $\mu$ L (0.02 ~ 0.1 mg/mL)