
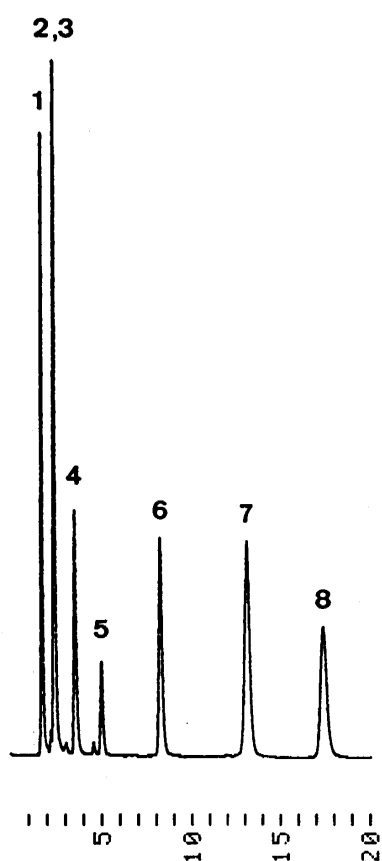
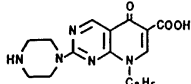
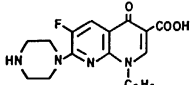
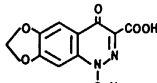
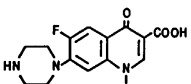
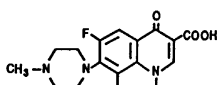
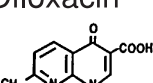
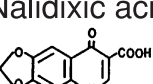


ナリジクス酸系抗生物質  
Nalidixic acid antibiotics

O950428D

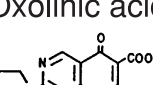

- 1   
Pipemidic acid
- 2   
Enoxacin
- 3   
Cinoxacin
- 4   
Norfloxacin
- 5   
Ofloxacin
- 6   
Nalidixic acid
- 7   
Oxolinic acid
- 8   
Piromidic acid

Column : YMC-Pack ODS-AM (5  $\mu$  m, 12nm)  
150  $\times$  4.6 mm I.D.

Eluent : 50mM  $\text{NH}_4\text{H}_2\text{PO}_4$  -  $(\text{NH}_4)_2\text{HPO}_4$  (pH 6.7)/tetrahydrofuran (90/10)  
containing 20mM TBA·Br

Flow rate : 1.0 mL/min

Temperature : ambient (25  $^\circ\text{C}$ )

Detection : UV at 260 nm, 0.32 AUFS

Injection : 5  $\mu$  L (0.03 mg/mL)