

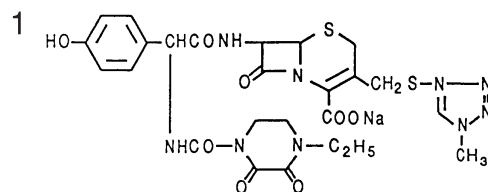


セフェム系抗生物質
Cephalosporin antibiotics


O940928S

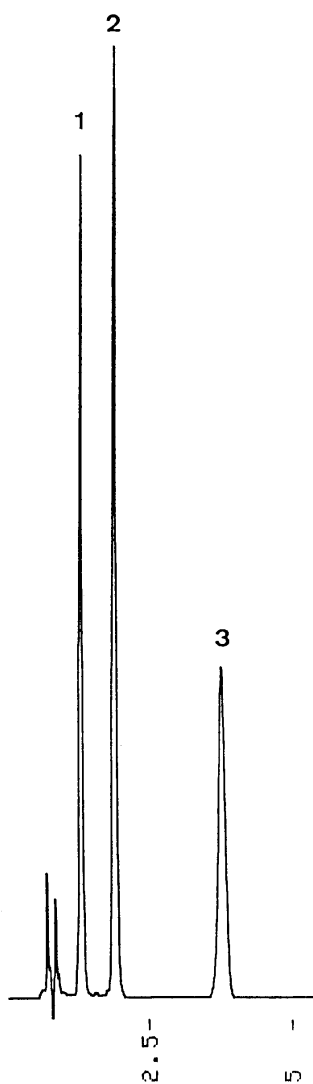

Cefoperazone sodium

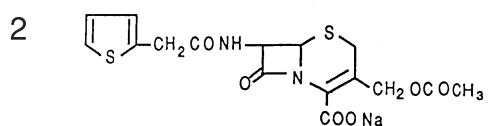

Cephalothin sodium

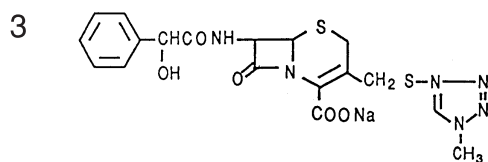

Cefamandole sodium

Column	: YMC-Pack ODS-A (3 μm, 12 nm) 50 X 4.6 mmI.D.
Eluent	: THF/50 mM NH ₄ H ₂ PO ₄ (15/85)
Flow rate	: 1.0 mL/min
Temperature	: 37°C
Detection	: UV at 260 nm, 0.16 AUFS
Injection	: 6 μL (0.17 mg/mL)