


セフェム系抗生物質
Cephalosporin antibiotics

O940719A

1. Cefsulodin sodium
2. Cephalexin hydrate
3. Cefaclor
4. Cephradine
5. Cephaloglycin
6. Cephapirin sodium
7. Cefotaxime sodium
8. Cephaloridine
9. Cefazolin sodium
10. Cefuroxime sodium
11. Cefoxitin sodium

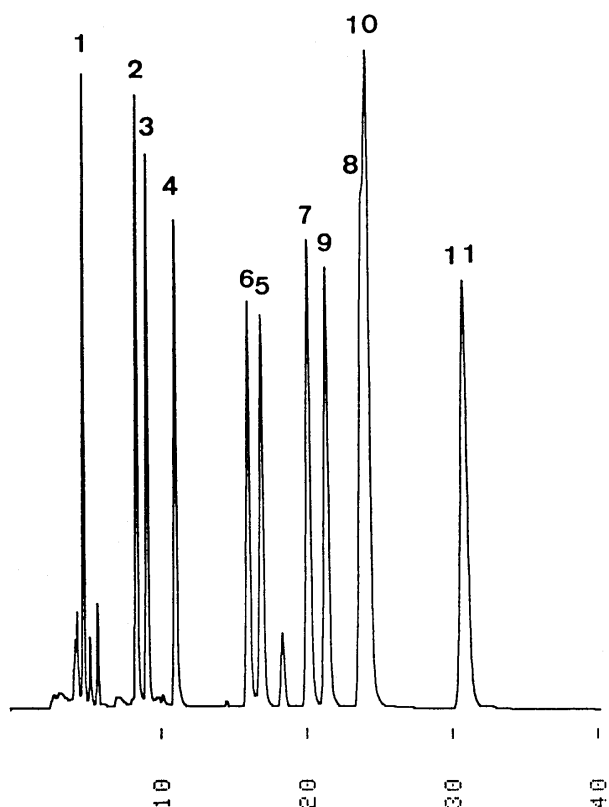

Column : J'sphere ODS-H80 (4 μ m, 8nm)
250 \times 4.6 mm I.D.
Eluent : acetonitrile/100mM KH₂PO₄ (10/90)
Flow rate : 0.8 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.16 AUFS
Injection : 5 μ L (0.12~0.35 mg/mL)