

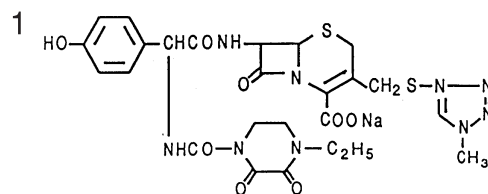


セフェム系抗生物質
Cephalosporin antibiotics


O940929D

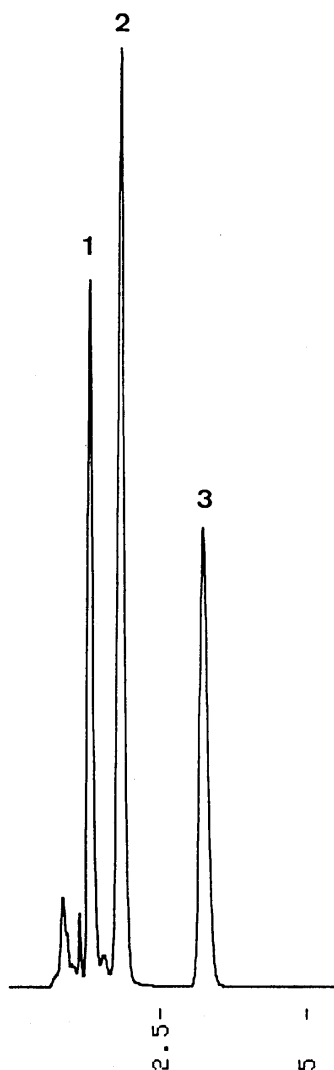

Cefoperazone sodium

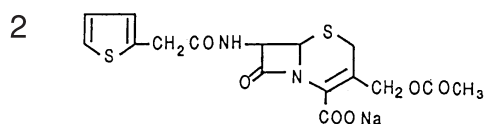

Cephalothin sodium

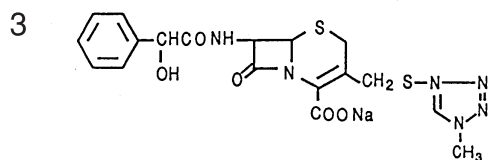

Cefamandole sodium

Column : J'sphere ODS-H80 (4 μ m, 8nm)
75 \times 2.0 mm I.D.
Eluent : THF/50mM NH₄H₂PO₄ (15/85)
Flow rate : 0.2 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.16 AUFS
Injection : 3 μ L (0.17 mg/mL)