


酸化防止剤 Antioxidants

A970526C

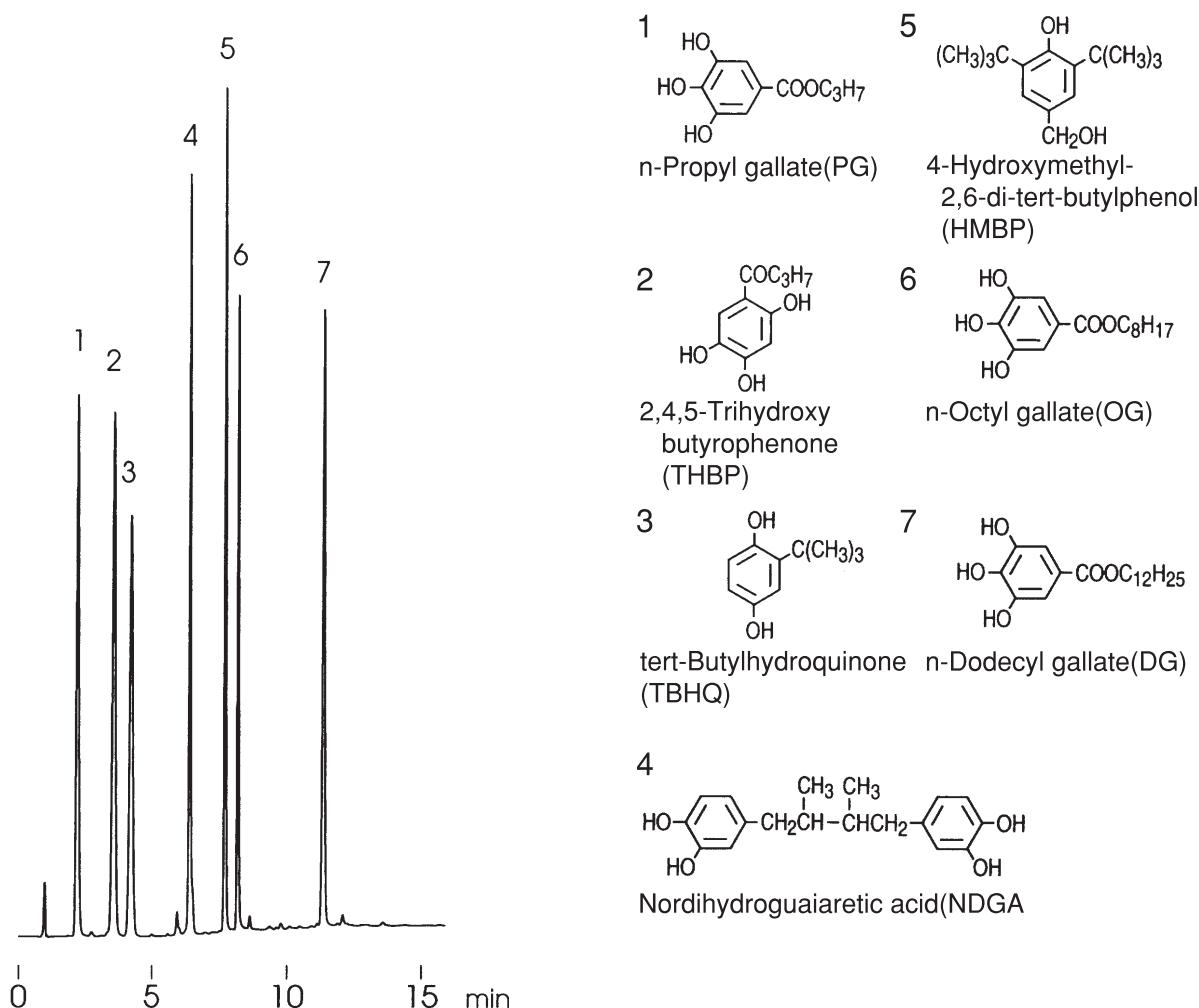

Column : YMC-Pack Pro C18 (5 μ m, 12nm)
 75 X 4.6 mm I.D.

Eluent : A)acetonitrile/methanol/5% acetic acid (20/20/60)
 B)acetonitrile/methanol/5% acetic acid (45/45/10)
 0-100% B (0-10min), 100% B (10-20min)

Flow rate : 1.0 mL/min

Temperature : 30 °C

Detection : UV at 280 nm, 0.32 AUFS

Injection : 9 μ L (0.05~0.6 mg/mL)