

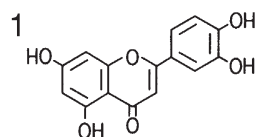


フラボノイド  
Flavonoids


S961107C

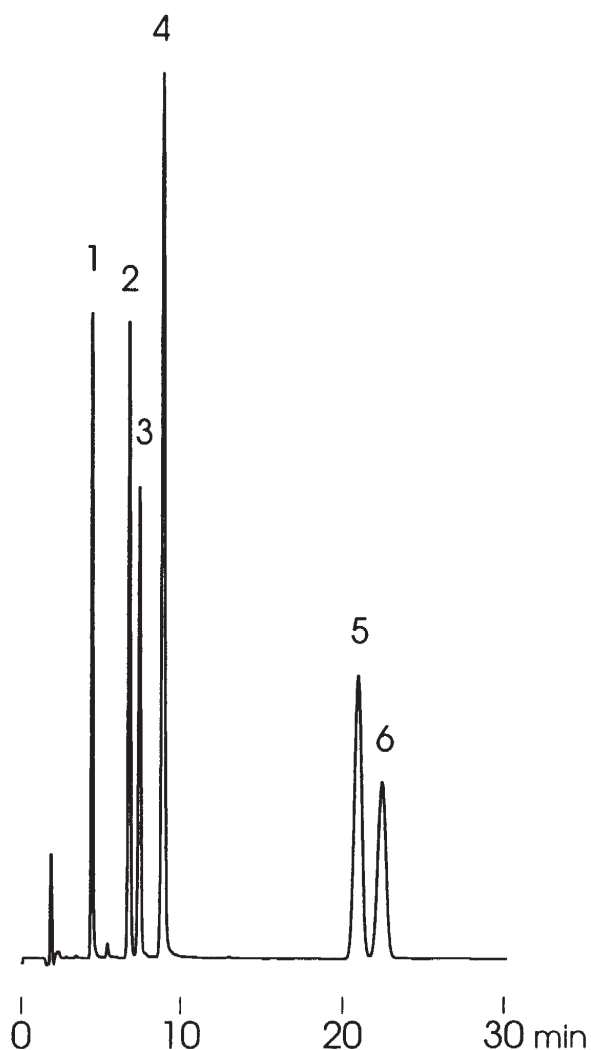


Luteolin

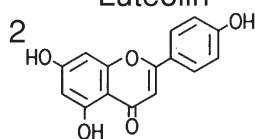


Apigenin

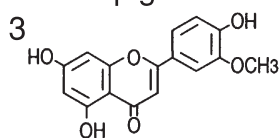


Chrysoeriol

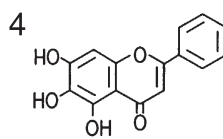

Baicalein

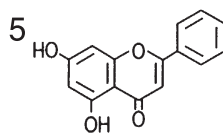

Chrysin

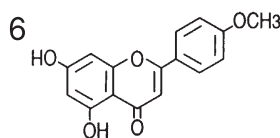

Acacetin

Column : YMC-Pack Pro C18 (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/10mM H<sub>3</sub>PO<sub>4</sub> (35/65)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 280 nm, 0.13 AUFS  
Injection : 5  $\mu$ L (0.05 mg/mL)