
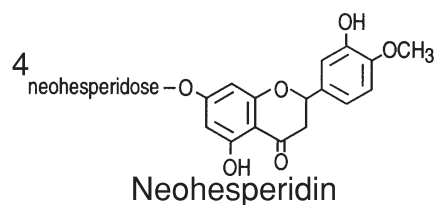
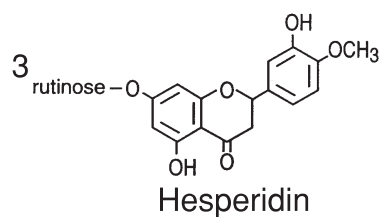
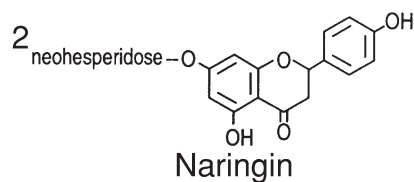
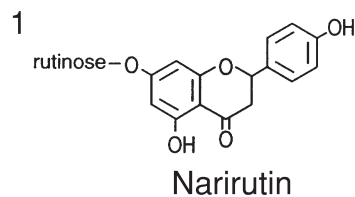


フラバノン配糖体  
Flavanone glycosides

A970506A

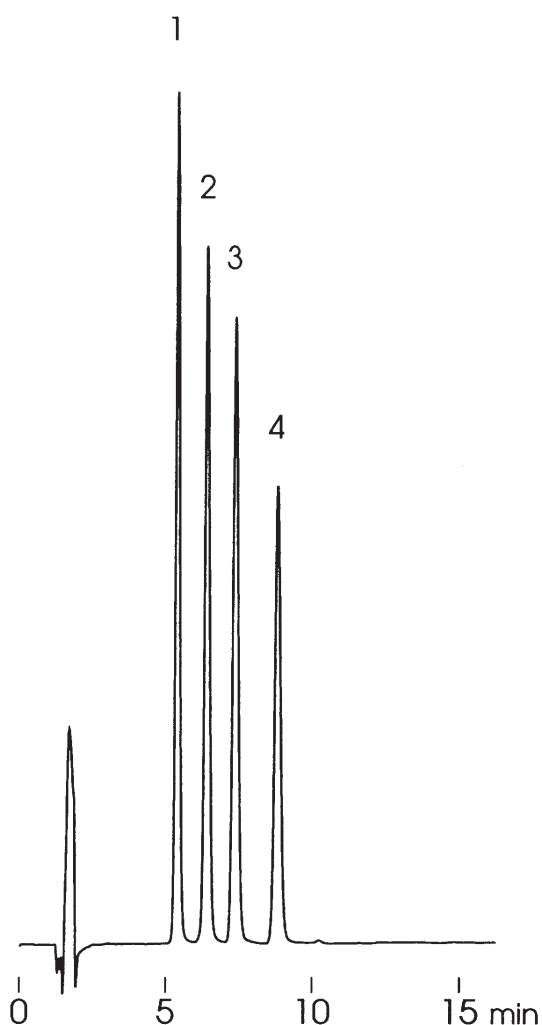

Column : J'sphere ODS-H80 (4  $\mu$  m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/10mM H<sub>3</sub>PO<sub>4</sub> (20/80)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 280 nm, 0.04 AUFS  
Injection : 7  $\mu$  L (0.025 mg/mL)