

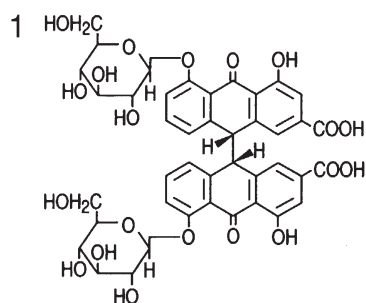


センノシド
Sennosides

A961126D

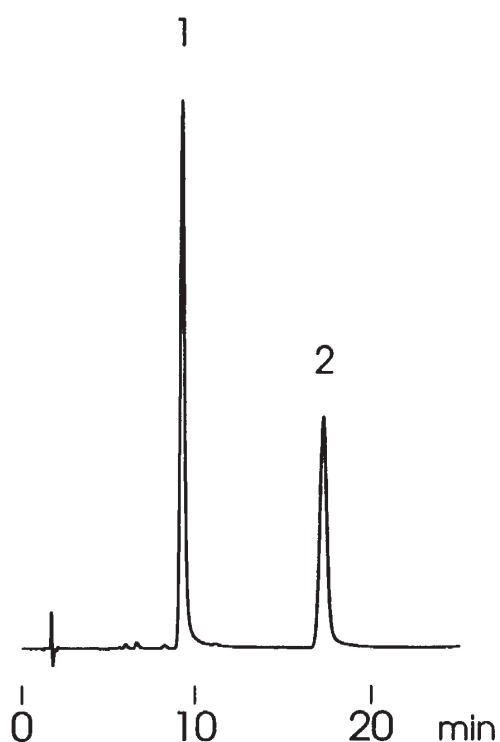

Sennoside B

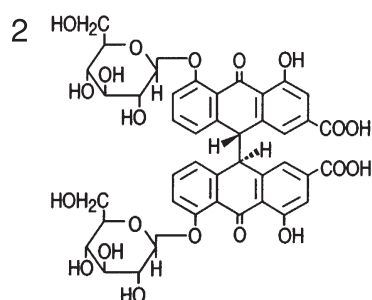

Sennoside A

| | |
|-------------|---|
| Column | : J'sphere ODS-M80 (4 μ m, 8nm) 150 \times 4.6 mm I.D. |
| Eluent | : methanol/water/acetic acid (30/65/5) |
| Flow rate | : 1.0 mL/min |
| Temperature | : 37 $^{\circ}$ C |
| Detection | : UV at 340 nm, 0.16 AUFS |
| Injection | : 10 μ L (0.1 mg/mL) |