
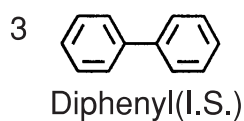
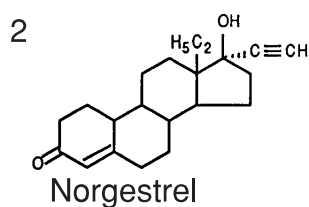
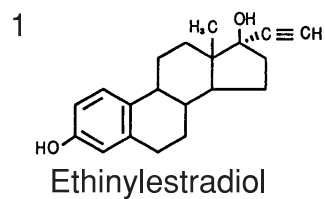


ノルゲストレルとエチニルエストラジオール  
Norgestrel and ethinylestradiol

S940427B

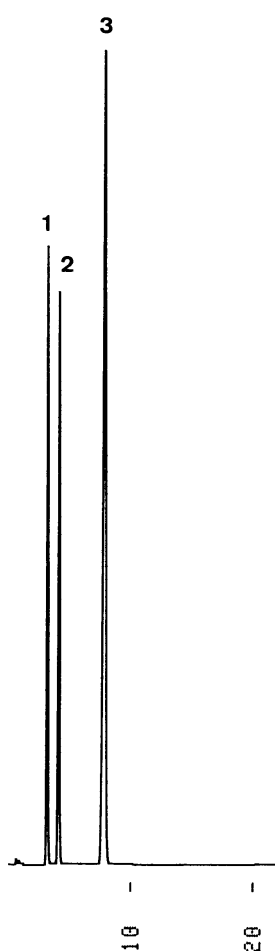

Column : J'sphere ODS-L80 (4  $\mu$  m, 8nm)  
75  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/water (50/50)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 230 nm, 0.13 AUFS  
Injection : 6  $\mu$  L (0.05~0.1 mg/mL)