


エクダイソン
Ecdysones

G910523H

1. β -Ecdysone
2. α -Ecdysone

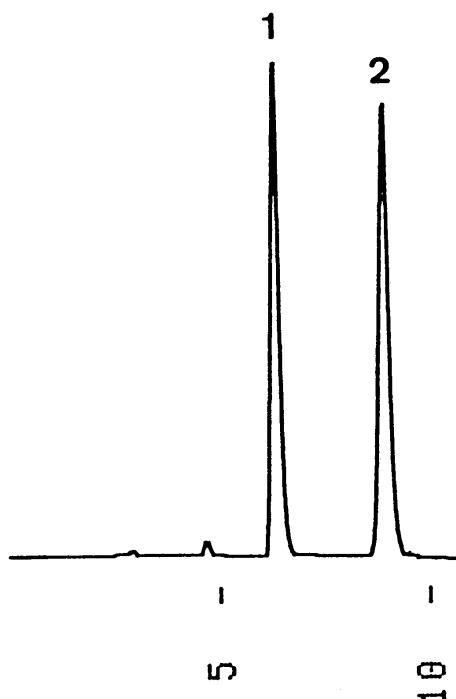

Courtesy of Dr.Furusawa, Kyoto Institute of Technology, Japan

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/methanol/water(5/50/45)
Flow rate : 0.5 mL/min
Temperature : ambient
Detection : UV at 254 nm, 0.04 AUFS
Injection : 5 μ L (0.04 mg/mL)