
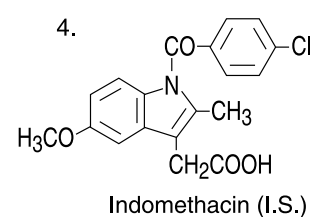
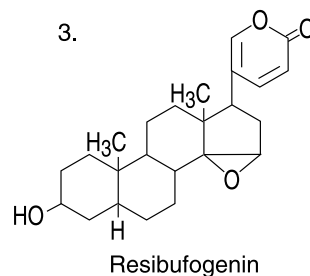
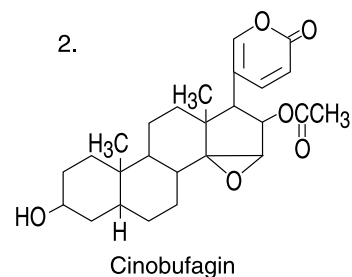
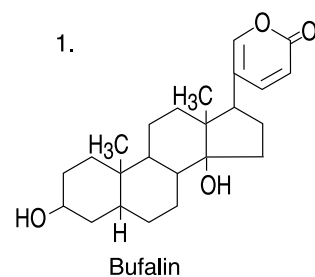


センソ成分

Steroids from toad venom

S010806A

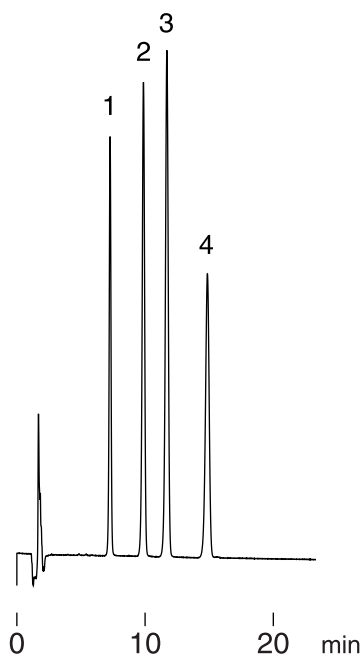

Column : YMC-Pack ODS-AL (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.

Eluent : acetonitrile/0.1% phosphoric acid (45/55)

Flow rate : 1.0 mL/min

Temperature : 40  $^{\circ}$ C

Detection : UV at 300 nm, 0.08 AUFS

Injection : 10  $\mu$ L (0.04~0.08 mg/mL)

(日本薬局方第14改正記載条件)