



鼻炎用点鼻剂  
Collunarium


Y920610A

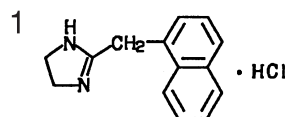

Naphazoline hydrochloride

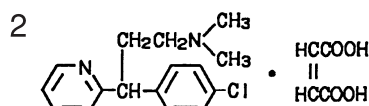

Chlorpheniramine maleate

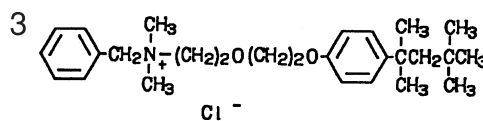

Benzethonium chloride

Courtesy of NIHON ARUTU PHARMACEUTICAL CO.,LTD.

Column	: YMC-Pack Ph (5 $\mu$ m, 12nm) 150 $\times$ 4.6 mm I.D.
Eluent	: acetonitrile/water/NH <sub>4</sub> SCN (48/52/0.5, v/v/w)
Flow rate	: 1.0 mL/min
Temperature	: 30 $^{\circ}$ C
Detection	: UV at 270 nm, 2.56 AUFS
Injection	: 10 $\mu$ L (0.2~5 mg/mL)
Sample	: collunarium