


塩化ベンゼトニウム  
Benzethonium chloride

T920527A

1. Benzethonium chloride

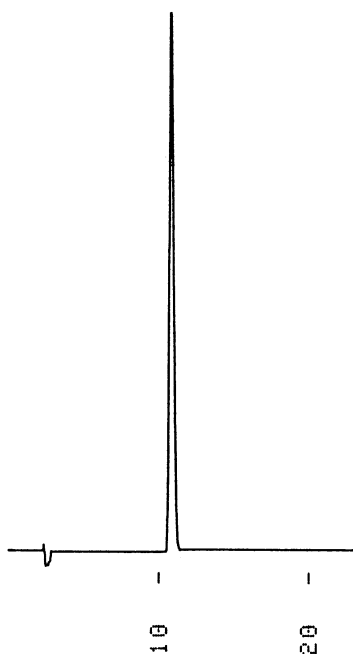

Courtesy of NIHON ARUTU PHARMACEUTICAL CO.,LTD.

Column	: YMC-Pack ODS-AM (5 $\mu$ m, 12nm) 150 $\times$ 6.0 mm I.D.
Eluent	: acetonitrile/water/NH <sub>4</sub> SCN (70/30/0.5, v/v/w)
Flow rate	: 1.0 mL/min
Temperature	: 30 $^{\circ}$ C
Detection	: UV at 270 nm, 0.16 AUFS
Injection	: 10 $\mu$ L (1.0 mg/mL)