
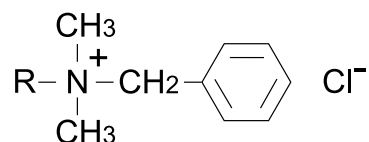


塩化ベンザルコニウム
Benzalkonium chloride

F031210G

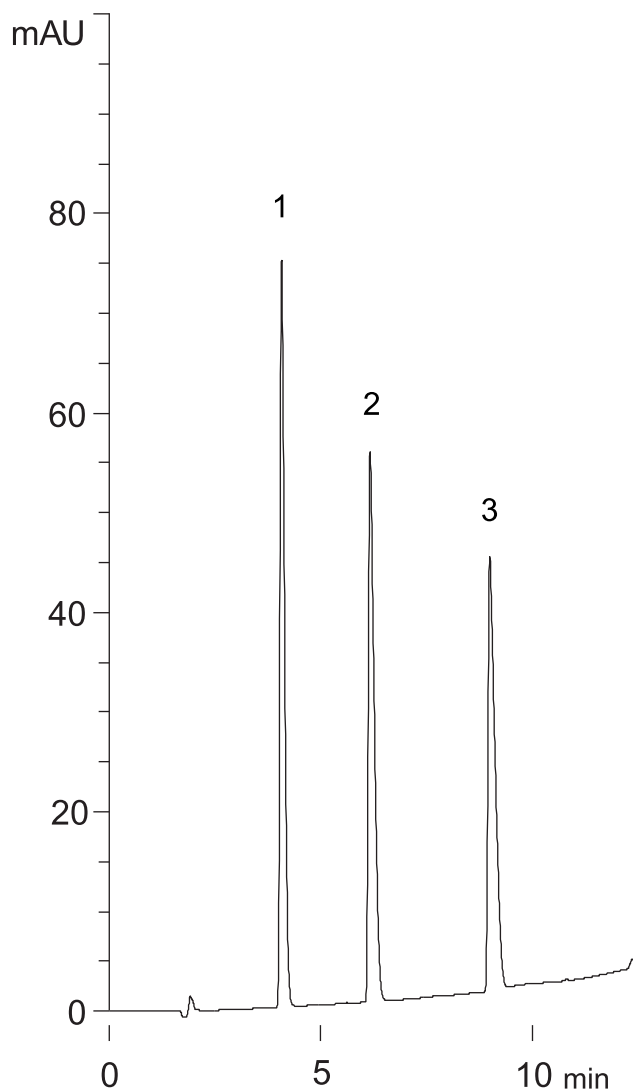

Benzalkonium chloride

1. R=-C₁₂H₂₅

2. R=-C₁₄H₂₉

3. R=-C₁₆H₃₃

Column : YMC-Pack Pro C18 (5 μm, 12nm)
150 × 4.6 mm I.D.
Eluent : A) 100mM NaClO₄/methanol (20/80)
B) methanol
0-50%B(0-10min)
Flow rate : 1.0 mL/min
Temperature : 37 °C
Detection : UV at 260 nm
Injection : 10 μL (1.0 mg/mL)