
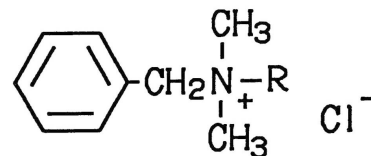


塩化ベンザルコニウム  
Benzalkonium chloride

T920529D

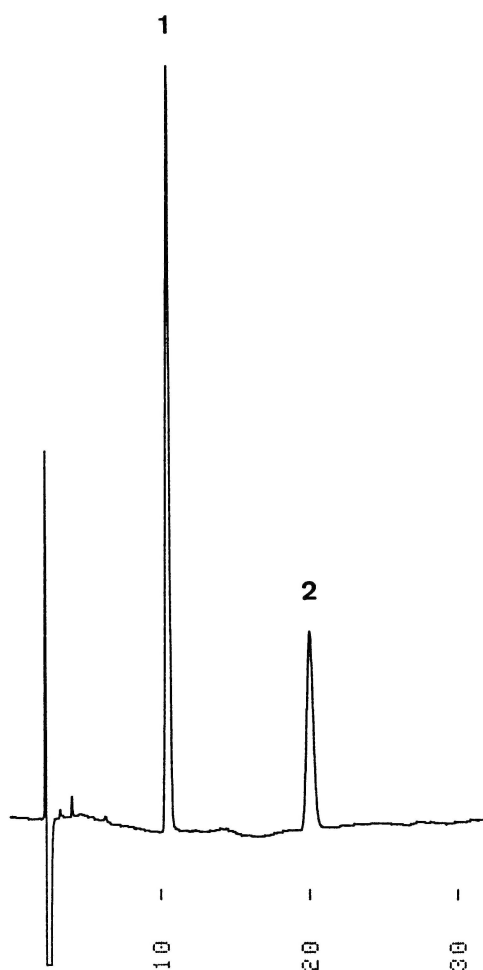

1. Benzalkonium chloride R=-C<sub>12</sub>H<sub>25</sub>
2. Benzalkonium chloride R=-C<sub>14</sub>H<sub>29</sub>

Courtesy of NIHON ARUTU PHARMACEUTICAL CO.,LTD.

Column	: YMC-Pack ODS-AM (5 μm, 12nm) 150 × 6.0 mm I.D.
Eluent	: acetonitrile/water/NH <sub>4</sub> SCN (70/30/0.5, v/v/w)
Flow rate	: 1.0 mL/min
Temperature	: 30 °C
Detection	: UV at 260 nm, 0.02 AUFS
Injection	: 10 μL (1.0 mg/mL benzalkonium chloride)