
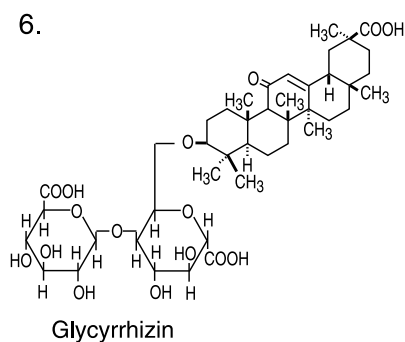
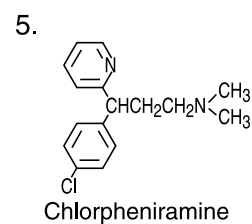
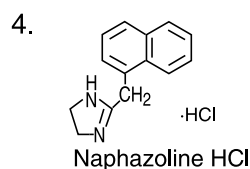
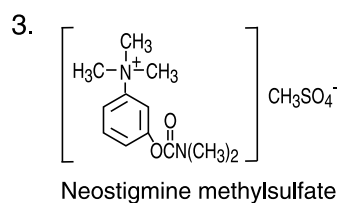
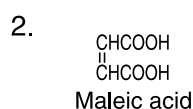
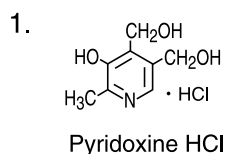


点眼剂成分
Ingredients in an eye drop

J010406A

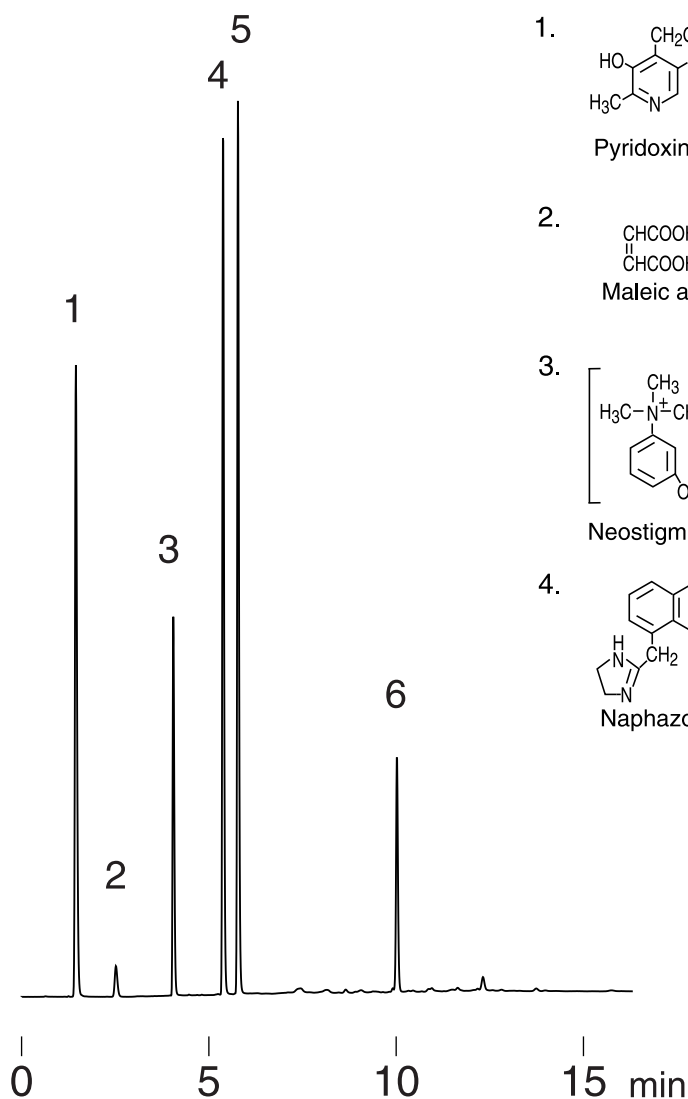

Column : Hydrosphere C18 (3 μ m, 12nm)
50 \times 4.6 mm I.D.
Eluent : A)20mM KH₂PO₄ - H₃PO₄ (pH 2.5)
B)acetonitrile
0-60%B(0-10min), 60%B(10-15min)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 265 nm, 0.128 AUFS
Injection : 5 μ L (0.03~0.3 mg/mL)