
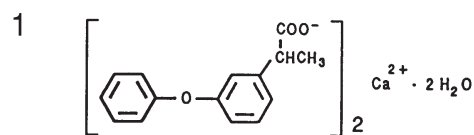



酸性非ステロイド性抗炎症薬

Acidic non-steroidal anti-inflammatory drugs


O950125C

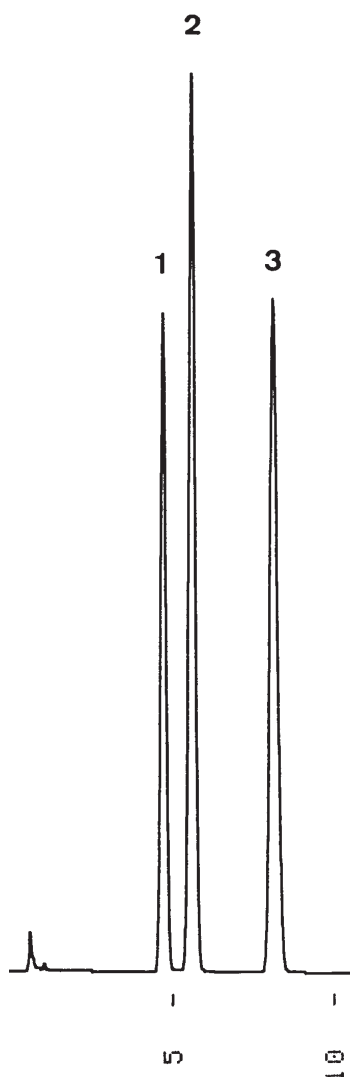

Fenopropfen calcium dihydrate

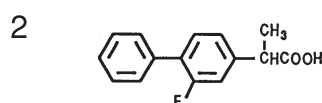

Flurbiprofen

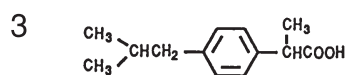

Ibuprofen

Column : J'sphere ODS-H80 (4 μ m, 8nm)
75 \times 4.6 mm I.D.
Eluent : methanol/water/acetic acid (65/35/1)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.04 AUFS
Injection : 7 μ L (160, 16, 820 μ g/mL)