
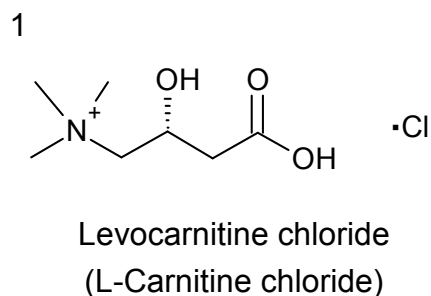




レボカルニチン注射液 (L-カルニチン注射液)
Levocarnitine injection (L-Carnitine injection)

V121023A

(A) Standard

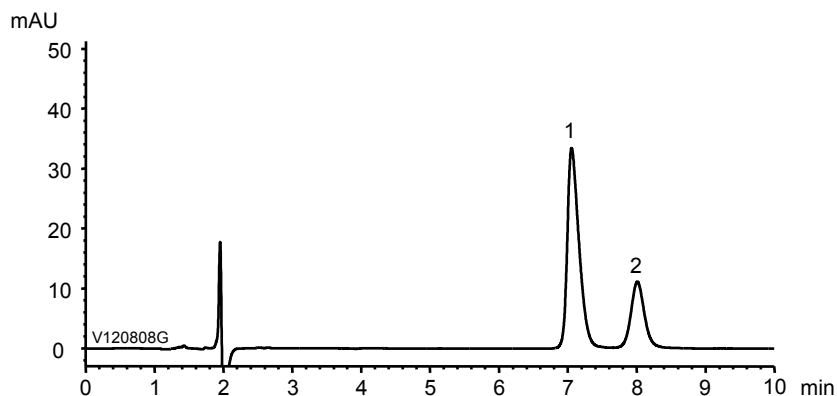

(B) Levocarnitine injection

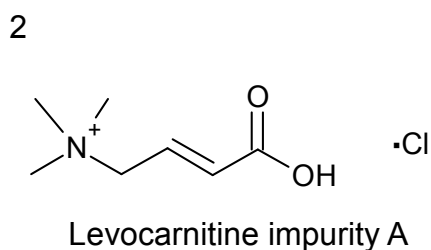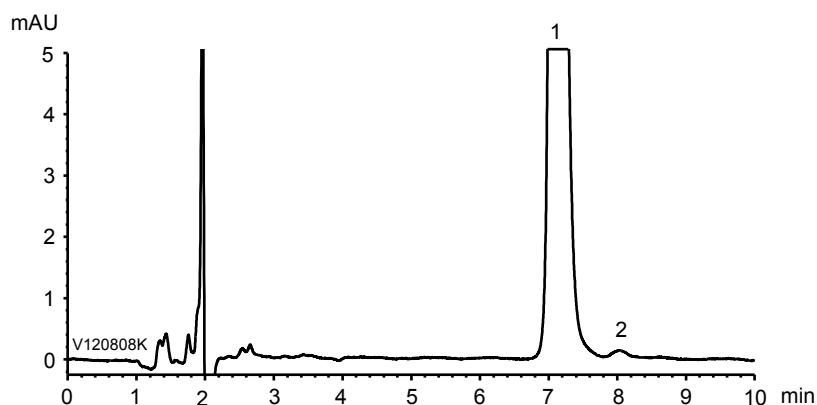

Column	: YMC-Triart Diol-HILIC (5 μm, 12 nm) 150 X 4.6 mmI.D.
Eluent	: 50 mM KH ₂ PO ₄ -H ₃ PO ₄ (pH 4.3)/acetonitrile (25/75)
Flow rate	: 1.0 mL/min
Temperature	: 30°C
Detection	: UV at 205 nm
Injection	: 2 μL
Sample	: (A) Standard (5 mg/mL, 0.025 mg/mL) (B) Levocarnitine injection diluted with acetonitrile/water (50/50) (correspond to 5 mg/mL Levocarnitine)