
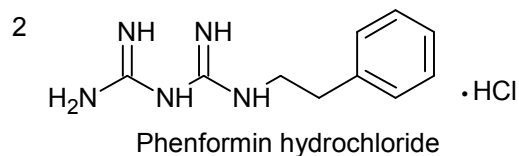
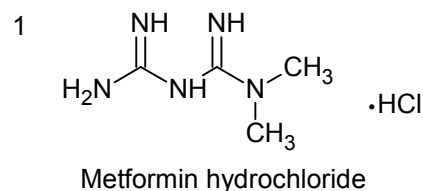


フェンホルミンおよびメトホルミン
Phenformin and metformin

U120420A

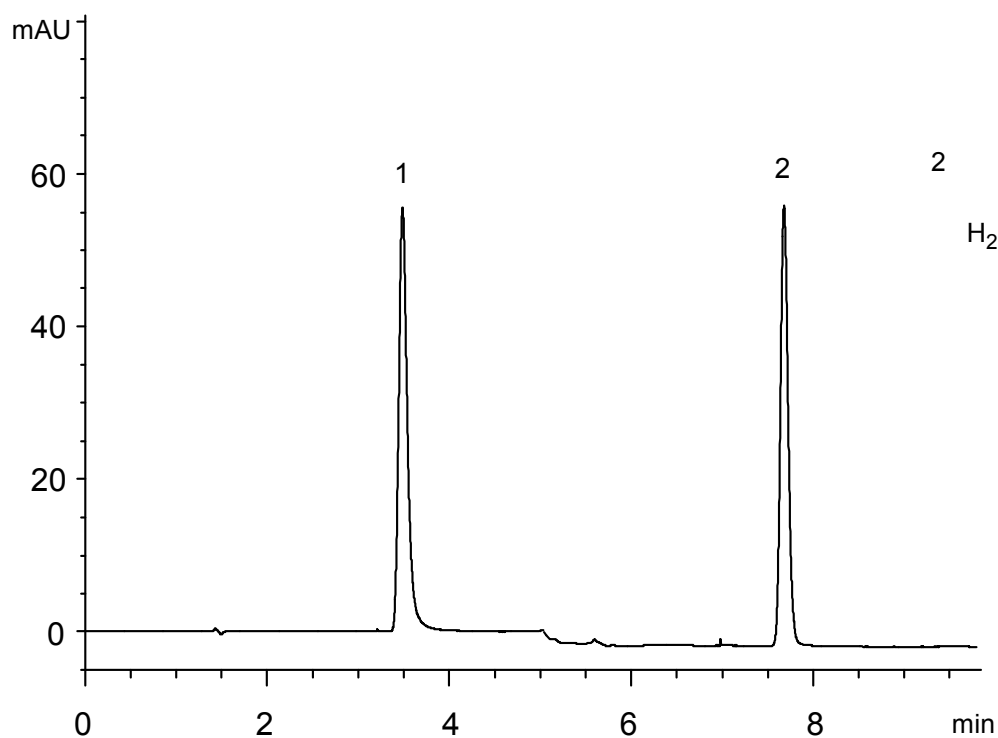

Column	: YMC-Triart C18 (5 μm 12 nm) 150 X 4.6 mm I.D.
Eluent	: A) 20 mM KH ₂ PO ₄ -K ₂ HPO ₄ (pH 7.0) B) acetonitrile 0%B(0-1 min), 0-20%B(1-3 min), 20%B(3-10 min)
Flow rate	: 1.0 mL/min
Temperature	: 25°C
Detection	: UV at 235 nm
Injection	: 10 μL (10 μg/mL)