

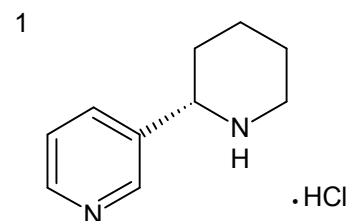


ニコチン及び関連化合物  
Nicotine and related compounds


U120409C

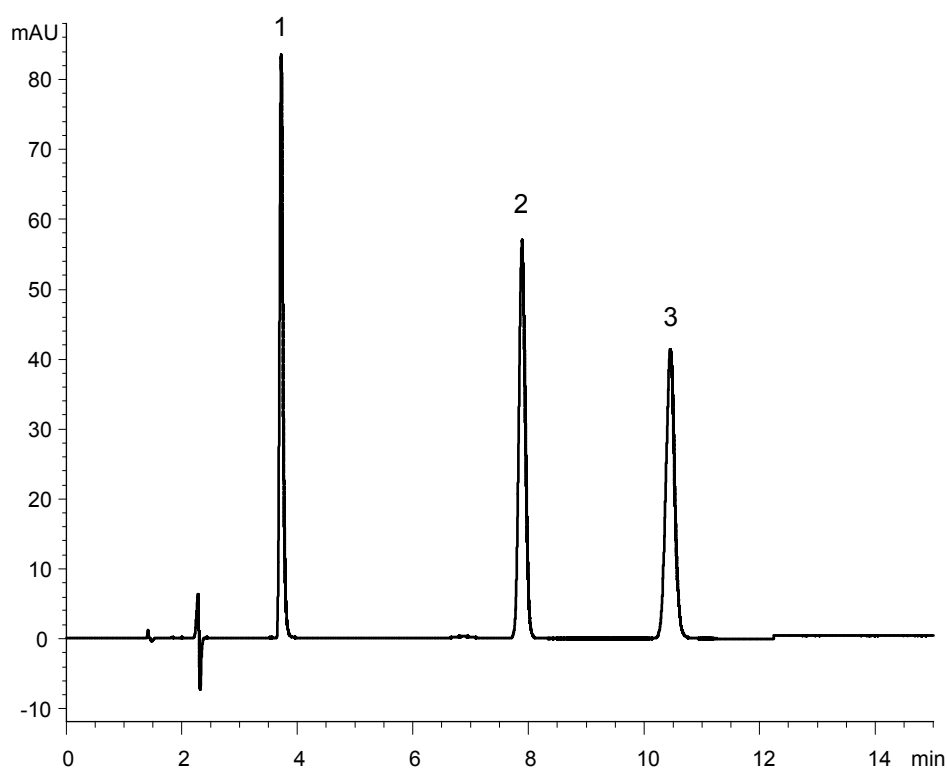

Anabasine hydrochloride

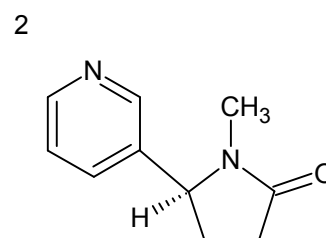

Cotinine

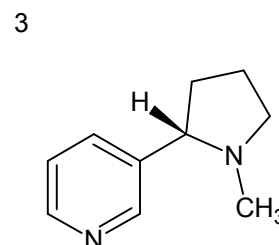

Nicotine

Column	: YMC-Triart C18 (3 μm, 12 nm) 150 X 3.0 mm I.D.
Eluent	: 20 mM KH <sub>2</sub> PO <sub>4</sub> -K <sub>2</sub> HPO <sub>4</sub> (pH 7.0)/acetonitrile (90/10)
Flow rate	: 0.425 mL/min
Temperature	: 40°C
Detection	: UV at 254 nm
Injection	: 2 μL (0.2 mg/mL)