



カルベジロール錠
Carvedilol Tablets


U110710A

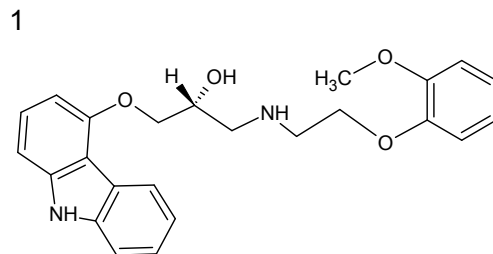

Carvedilol

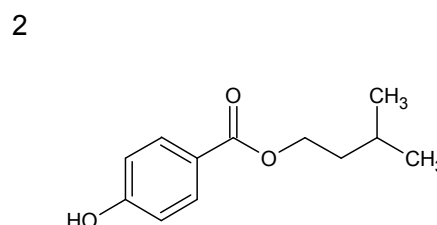

Isoamyl *p*-hydroxybenzoate (I.S.)

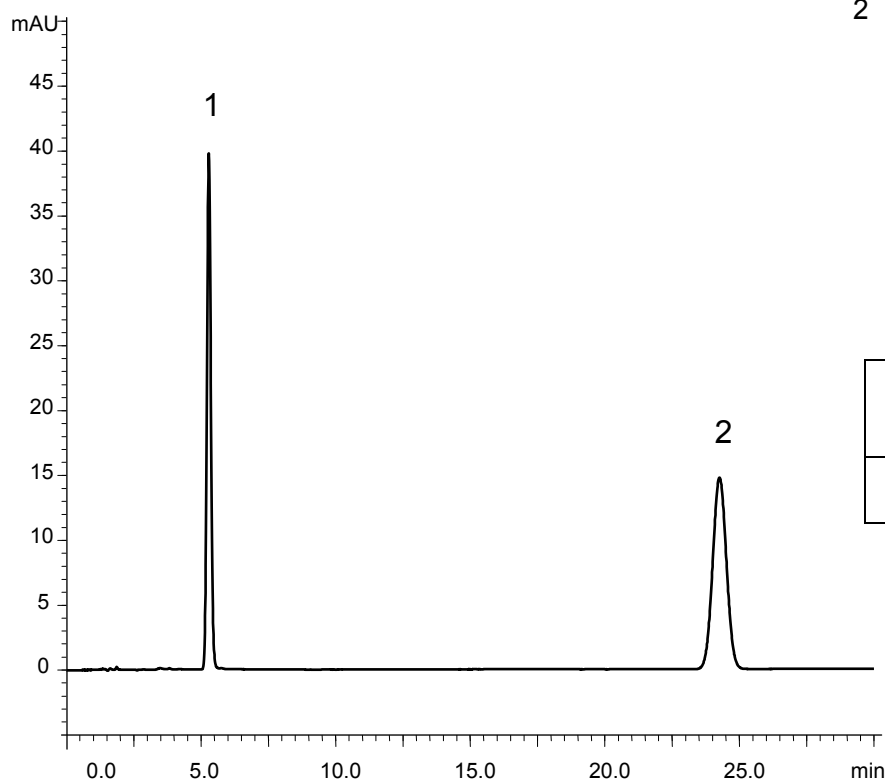

	System suitability requirement	result
Resolution (1, 2)	≥20	31.2

Column : YMC-Triart C18 (5 μm, 12 nm)
150 X 4.6 mm I.D.

Eluent : 20 mM KH₂PO₄-K₂HPO₄ (pH 5.0)/methanol (45/55)

Flow rate : 1.2 mL/min (*adjust the flow rate so that the retention time of carvedilol is about 5 min*)

Temperature : 40°C

Detection : UV at 240 nm

Injection : 10 μL (0.01, 0.029 mg/mL)

(The Japanese Pharmacopoeia 16th; Assay)