
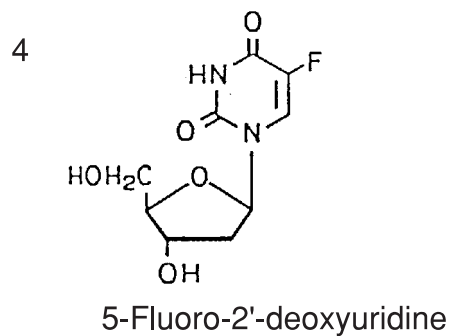
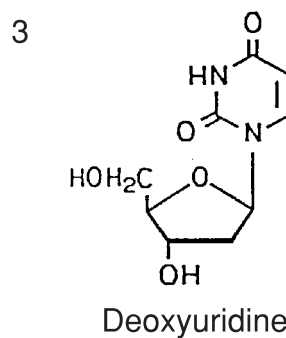
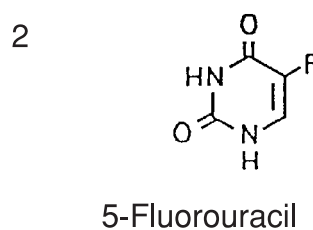
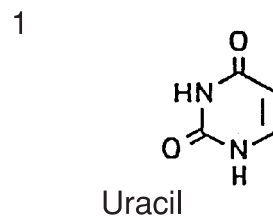


制がん剤および関連化合物

Antineoplastic and related compounds

T920924G

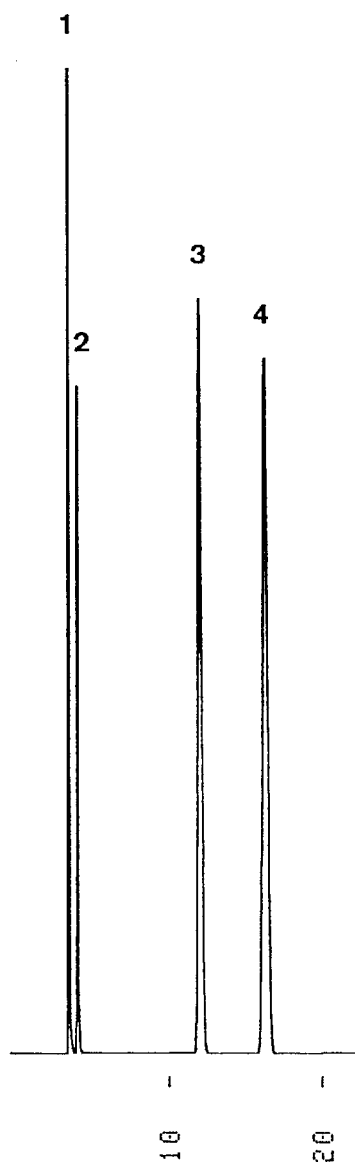

Column : YMC-Pack ODS-AQ (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 10mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub>  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 254 nm, 0.16 AUFS  
Injection : 10  $\mu$ L (0.025~0.2 mg/mL)