



制がん剤および関連化合物
Antineoplastic and related compounds

T920924E

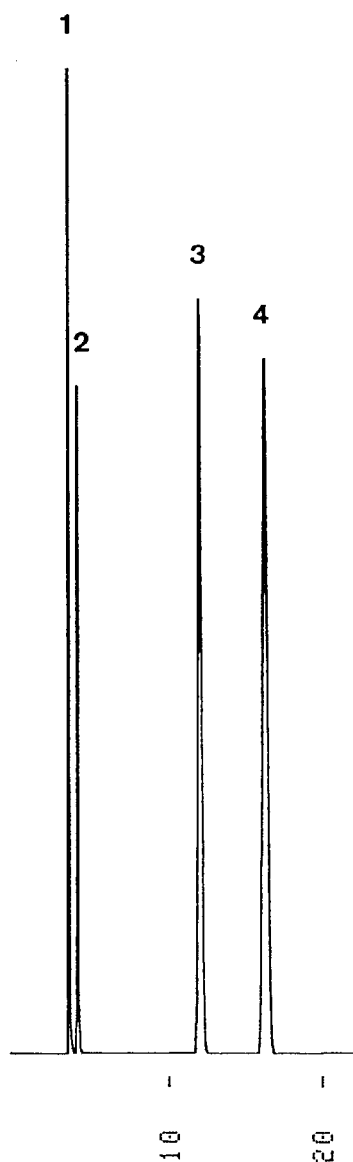


1.

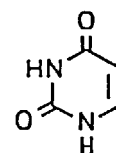


Uracil

2.

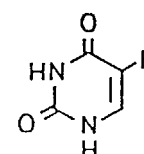


5-Fluorouracil

3.

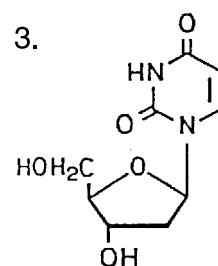

Deoxyuridine

4.

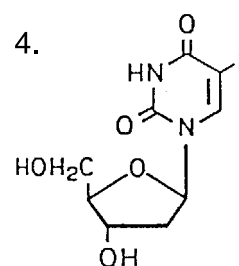

5-Fluoro-2'-deoxyuridine

Column : YMC-Pack ODS-A (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 10mM NH₄H₂PO₄
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 254 nm, 0.16 AUFS
Injection : 10 μ L (0.025~0.2 mg/mL)