

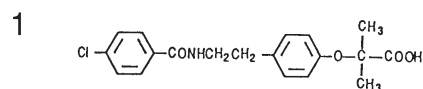


抗高血脂症剂  
Hypolipidemic drugs


S930706C

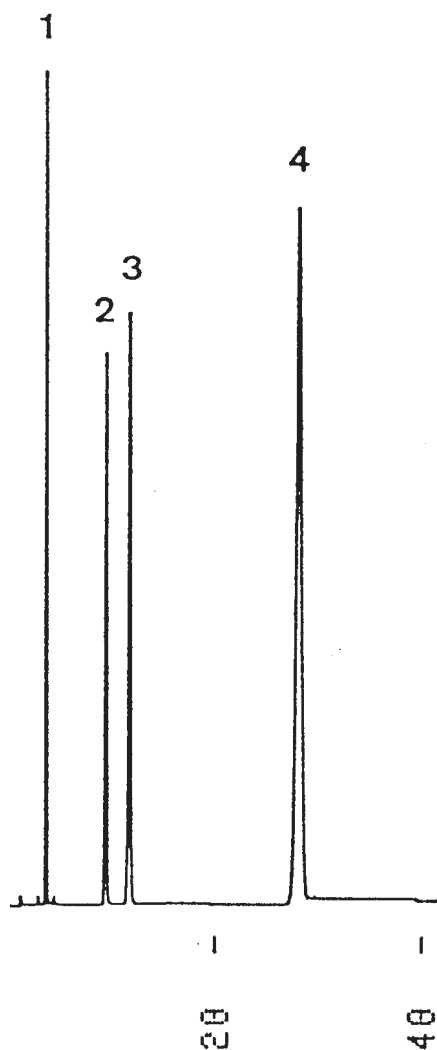


Bezafibrate

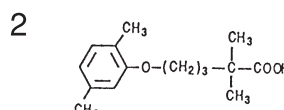

Gemfibrozil

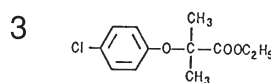

Clofibrate

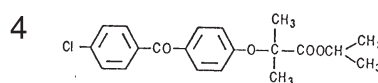

Fenofibrate

Column : J'sphere ODS-M80 (4  $\mu$  m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/water/acetic acid (60/40/1)  
Flow rate : 1.0 mL/min  
Temperature : 30  $^{\circ}$ C  
Detection : UV at 275 nm, 0.064 AUFS  
Injection : 8  $\mu$  L (0.05~0.25 mg/mL)