

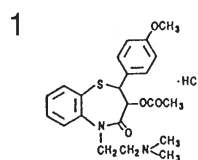


カルシウム拮抗薬
Calcium antagonists


S930702C

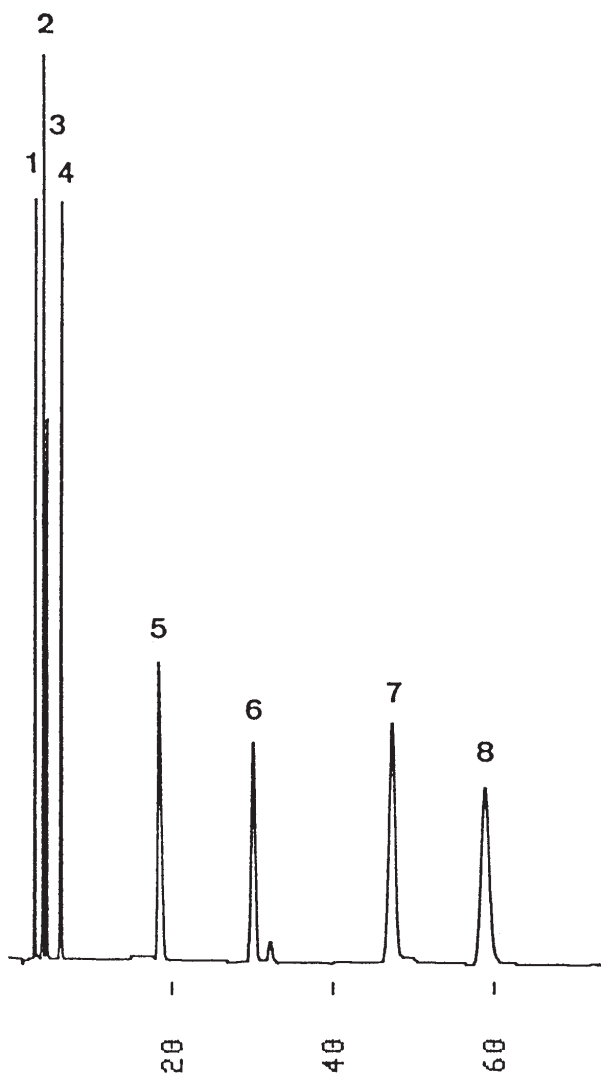


Diltiazem HCl

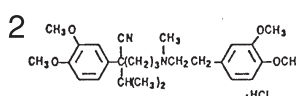


Verapamil HCl

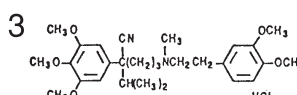


Methoxyverapamil HCl

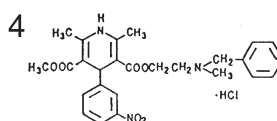


Nicardipine HCl

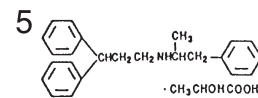


Prenylamine lactate

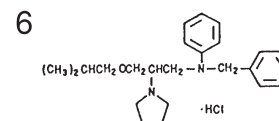

Bepridil HCl

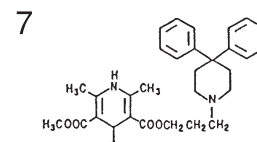

Nifedipine HCl

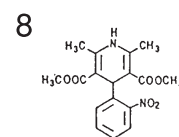

Nifedipine

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : 20mM KH₂PO₄ - H₃PO₄ (pH 3.0)/acetonitrile/THF/TBA·Br
(75/15/10/0.5, v/v/v/w)
Flow rate : 0.8 mL/min
Temperature : 30 °C
Detection : UV at 270 nm, 0.064 AUFS
Injection : 8 μ L (0.05~1.4 mg/mL)