

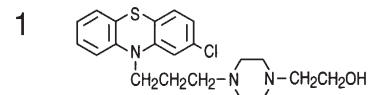


抗精神病薬  
Antipsychotic drugs


O950627B

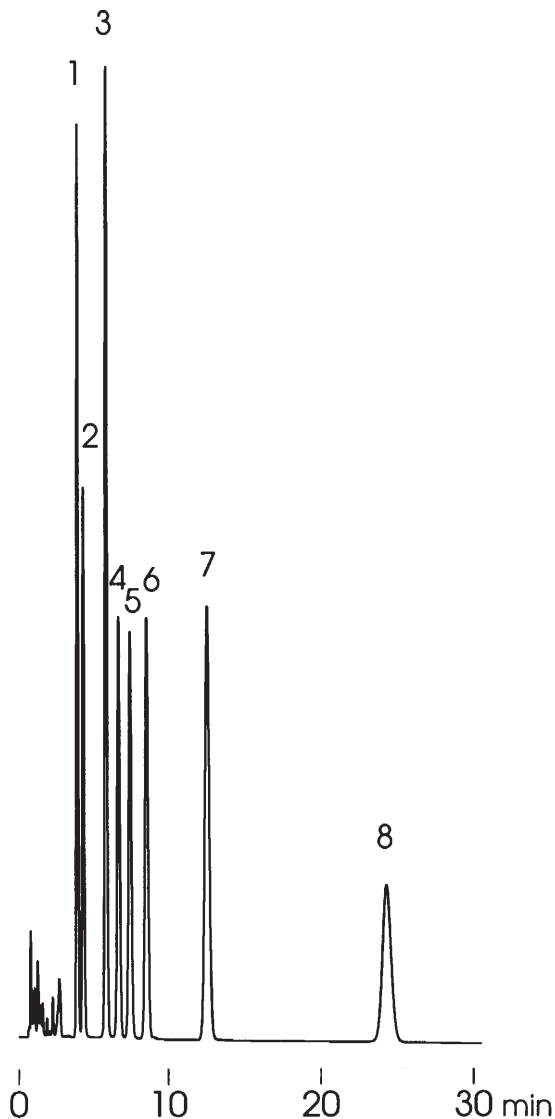


Perphenazine

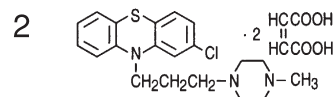


Prochlorperazine maleate

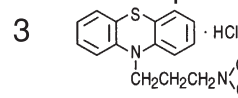


Promazine hydrochloride

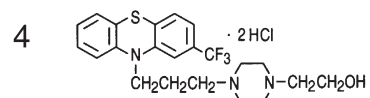


Fluphenazine hydrochloride

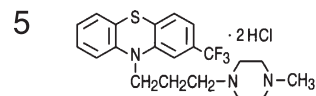


Trifluoperazine hydrochloride

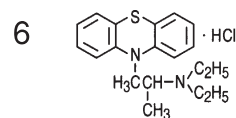

Ethopropazine hydrochloride  
(Profenamine hydrochloride)

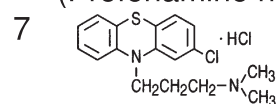

Chlorpromazine hydrochloride

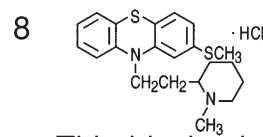

Thioridazine hydrochloride

Column : J'sphere ODS-H80 (4  $\mu$  m, 8nm)  
75  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/water/TFA (30/70/0.1)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 254 nm, 0.064 AUFS  
Injection : 5  $\mu$  L (0.02 mg/mL)

(Profenamine hydrochloride)