
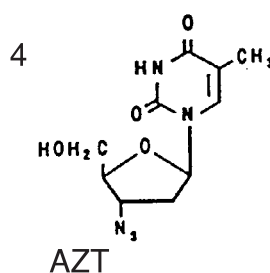
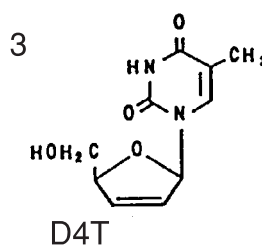
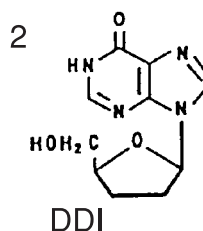
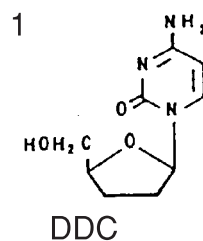


抗HIVヌクレオシド誘導体

Anti-human immunodeficiency virus (HIV) agents

K930311A

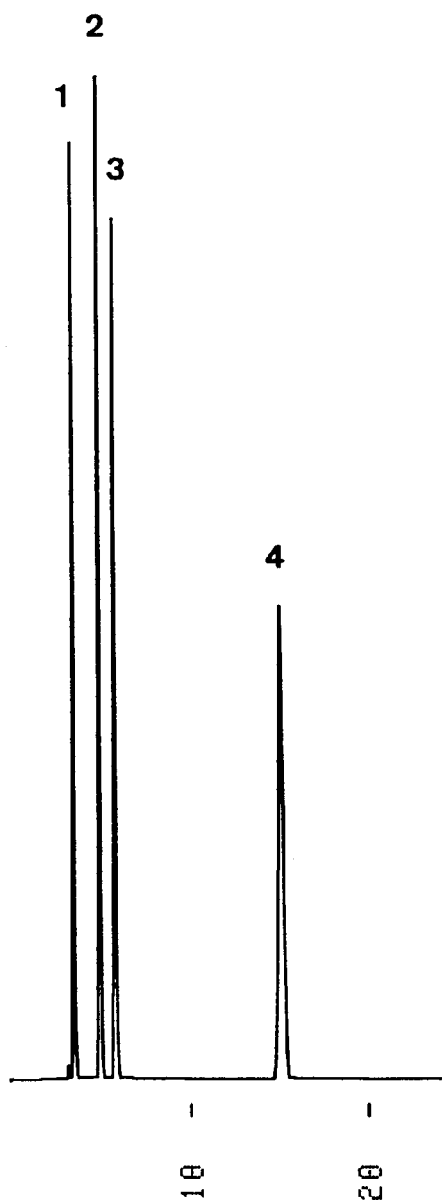

Column : YMC-Pack C8 (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : methanol/10mM KH<sub>2</sub>PO<sub>4</sub> (10/60)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 254 nm, 0.16 AUFS  
Injection : 7  $\mu$ L (0.125 mg/mL)