
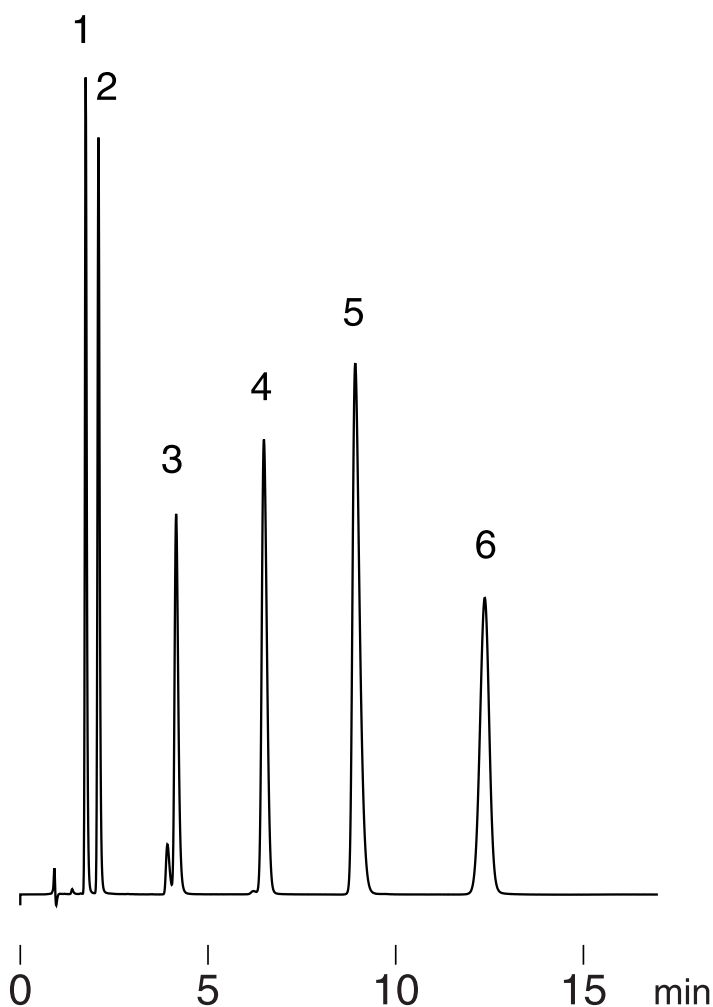
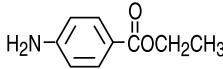
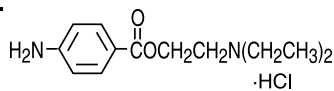
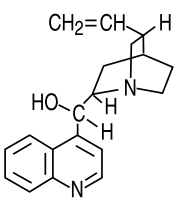
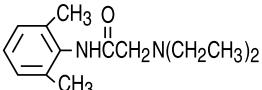
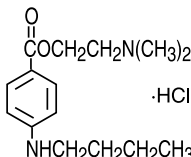


局所麻醉薬
Local anesthetics

J010302E

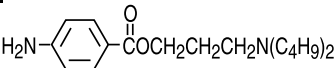

1. 
Benzocaine
(Ethyl 4-aminobenzoate)
2. 
Procaine HCl
3. 
Cinchonine
4. 
Lidocaine
5. 
Tetracaine HCl
6. 
Butacaine

Column : Hydrosphere C18 (3 μ m, 12nm)
75 \times 4.6 mm I.D.
Eluent : 20mM KH₂PO₄-K₂HPO₄ (pH 6.9) / methanol (40/60)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 254 nm, 0.16 AUFS
Injection : 10 μ L (0.02~0.5 mg/mL)