
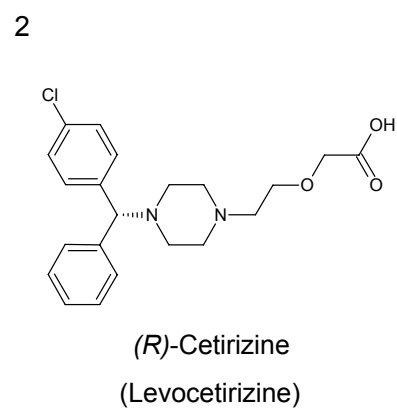
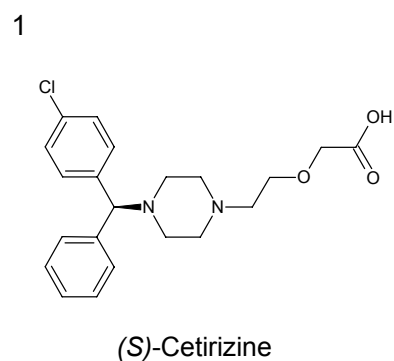


セチリジン (レボセチリジン)

Cetirizine (Levocetirizine)

F141211A

$\alpha$	1.3
$R_s$	3.5

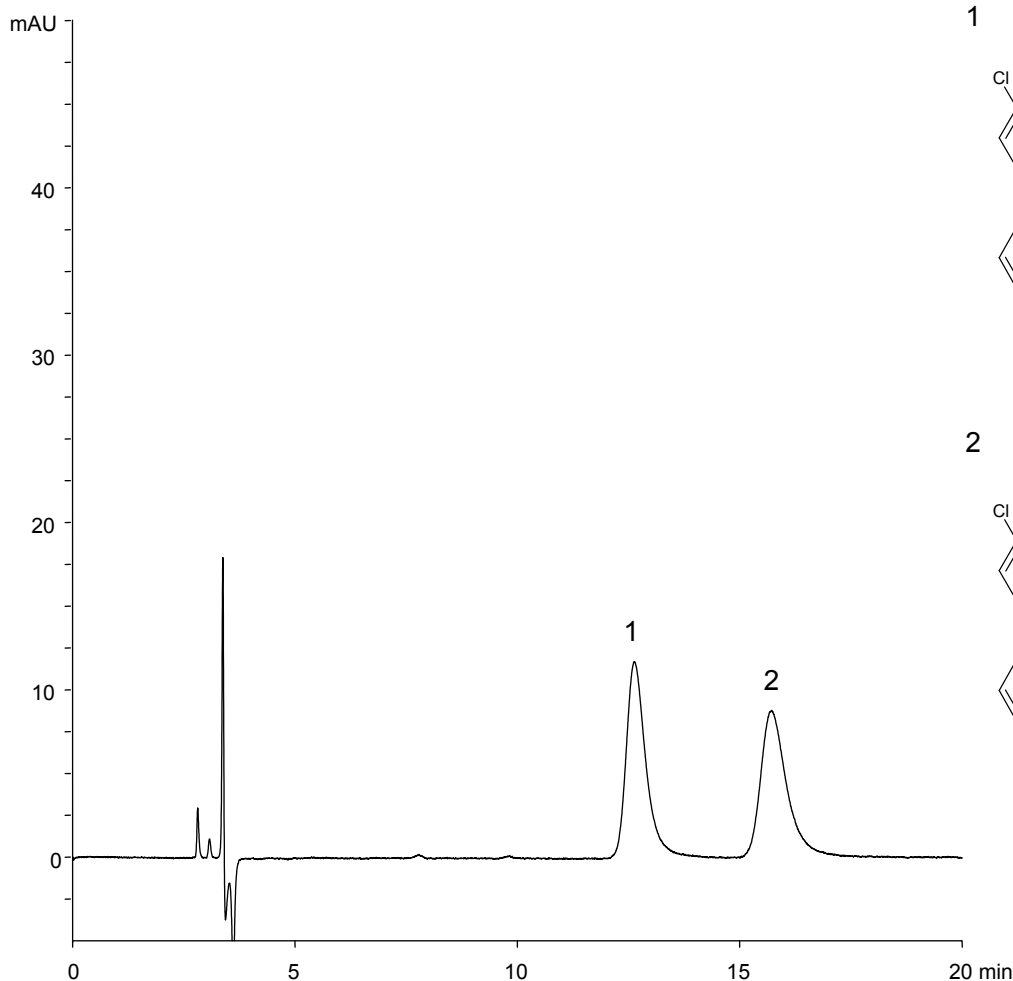

Column : CHIRAL ART Cellulose-C (5  $\mu$ m)  
250 X 4.6 mmI.D.

Eluent : *n*-hexane/2-propanol/formic acid/diethylamine (70/30/0.1/0.1)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 230 nm

Injection : 2  $\mu$ L (200  $\mu$ g/mL)