



デュロキセチン塩酸塩
Duloxetine hydrochloride


F131003A

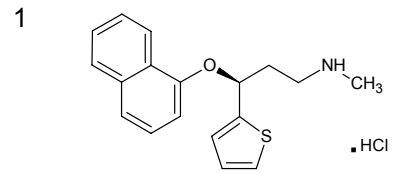

Duloxetine hydrochloride
(*S*)-Duloxetine hydrochloride

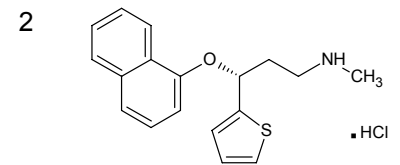

Duloxetine related compound A
(*R*)-Duloxetine hydrochloride

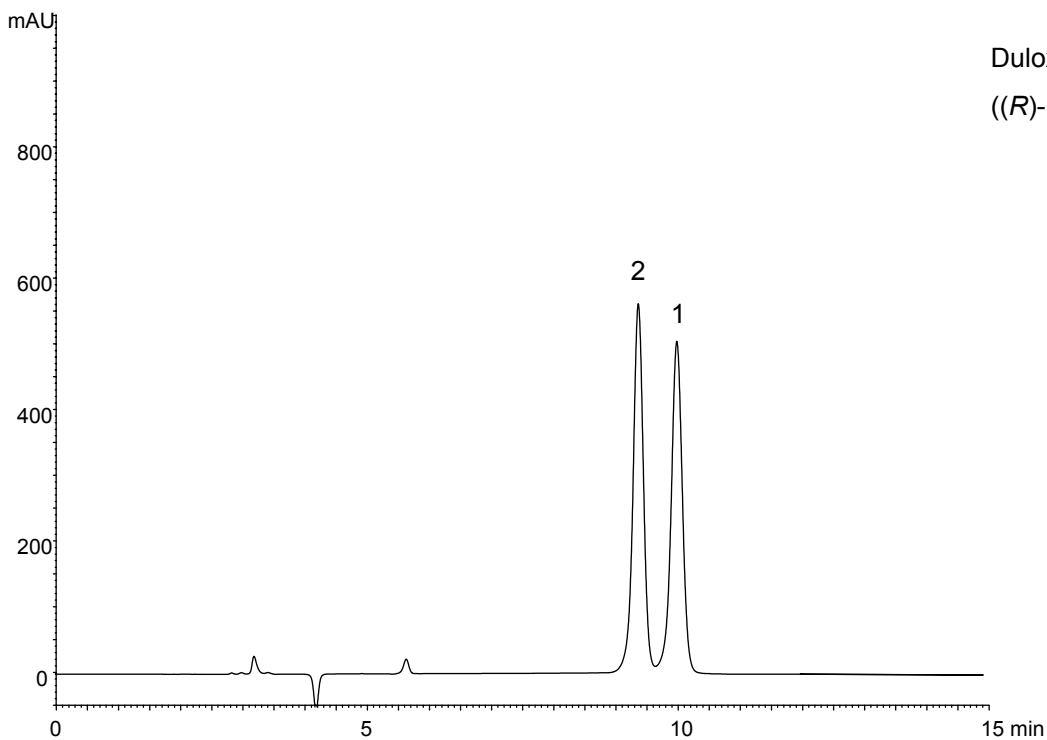

Column	: CHIRAL ART Amylose-C (5 μ m) 250 X 4.6 mm I.D.
Eluent	: <i>n</i> -hexane/ethanol/diethylamine (95/5/0.2)
Flow rate	: 1.0 mL/min
Temperature	: 30°C
Detection	: UV at 230 nm
Injection	: 10 μ L (0.1 mg/mL)