
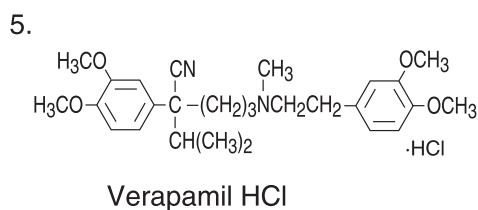
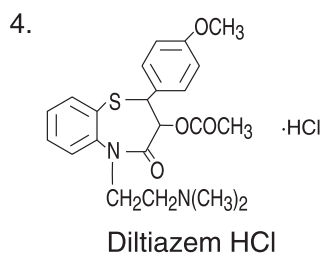
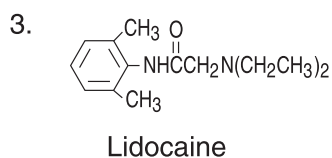
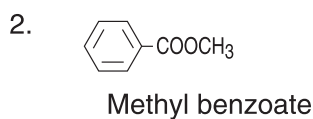
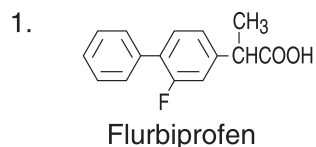


酸性および塩基性医薬品(pH 8.8)  
Acidic and basic drugs (pH 8.8)

F040714C

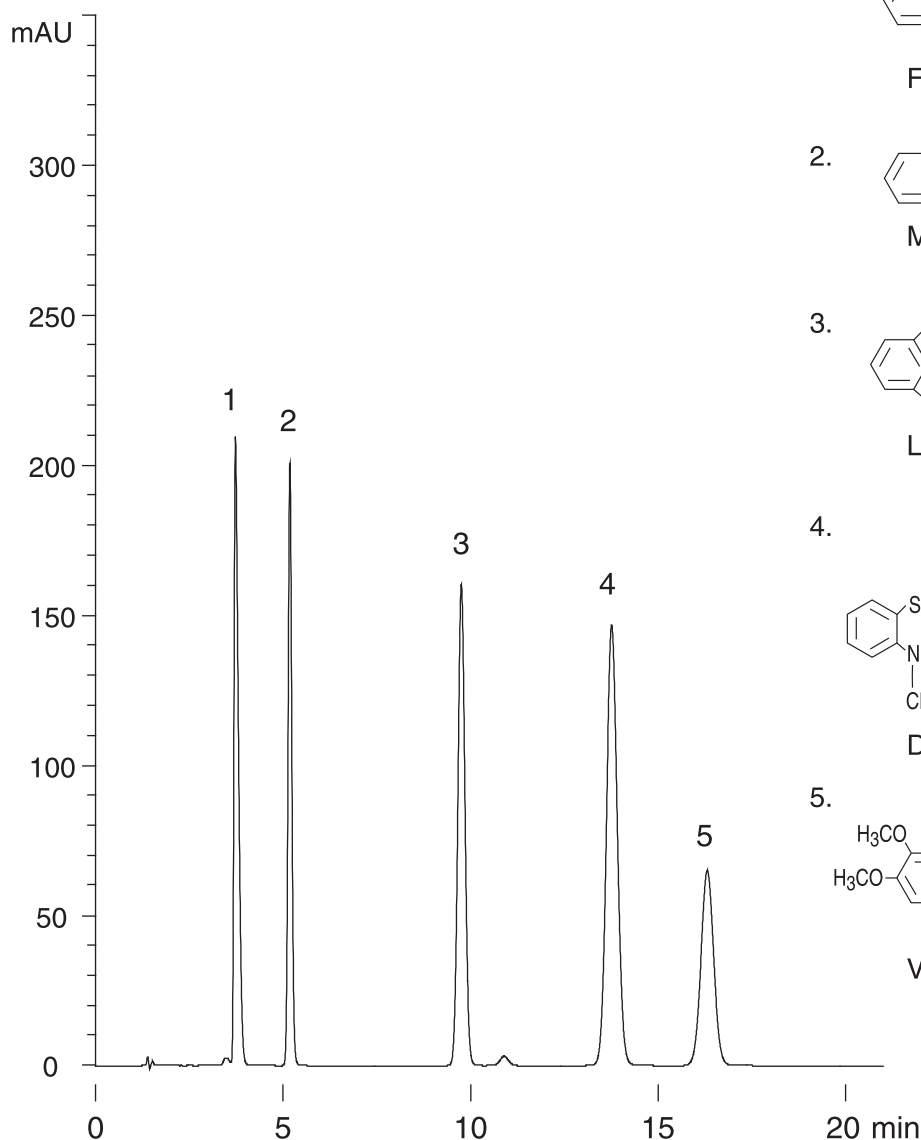

Column : YMC-Pack Pro C18 RS (5 μm, 8nm)  
150 × 4.6 mm I.D.  
Eluent : 20mM CH<sub>3</sub>COONH<sub>4</sub>-NH<sub>3</sub>(pH 8.8)/methanol (35/65)  
Flow rate : 1.0 mL/min  
Temperature : 30 °C  
Detection : UV at 220 nm  
Injection : 10 μL (0.04~100 μg/mL)