
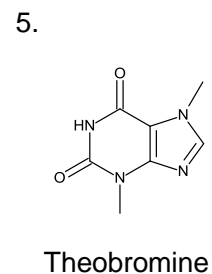
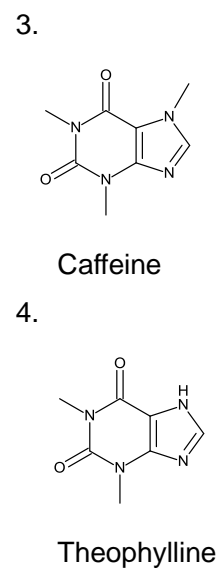
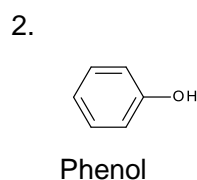
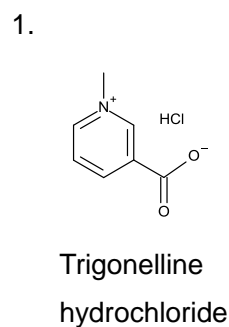


プリンアルカロイド  
Purine alkaloids

X150714A

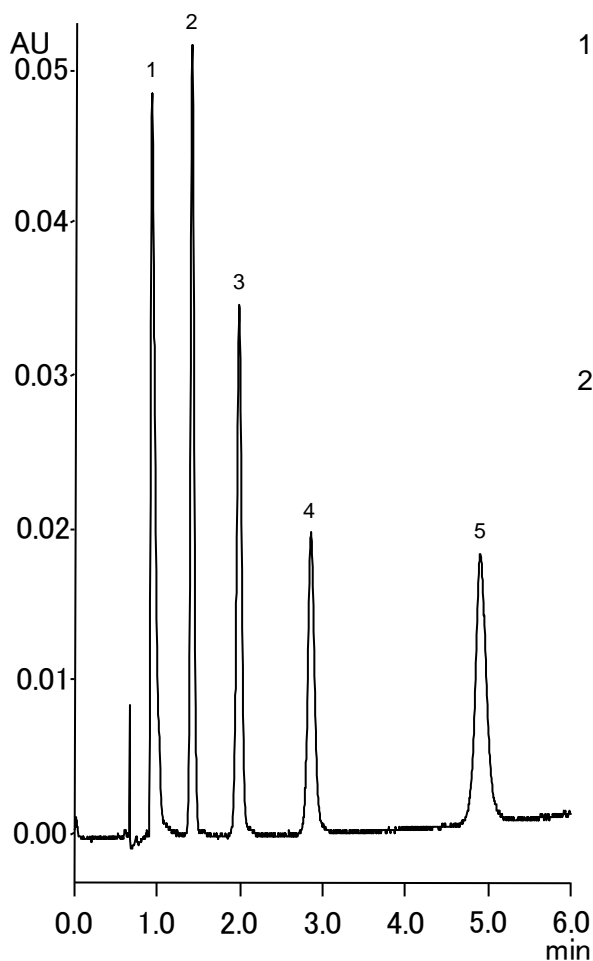

Column : Alcyon SFC Triart Diol (5  $\mu$ m, 12 nm)  
150 X 4.6 mm I.D.  
Eluent : CO<sub>2</sub>/ethanol (90/10)  
Flow rate : 3.0 mL/min  
Temperature : 40°C  
Detection : UV at 230 nm  
Back Pressure : 13.8 MPa (2000 psi)  
Injection : 5  $\mu$ L (0.085 mg/mL ~ 5.7 mg/mL)