



ステロイド
Steroids

X150801J

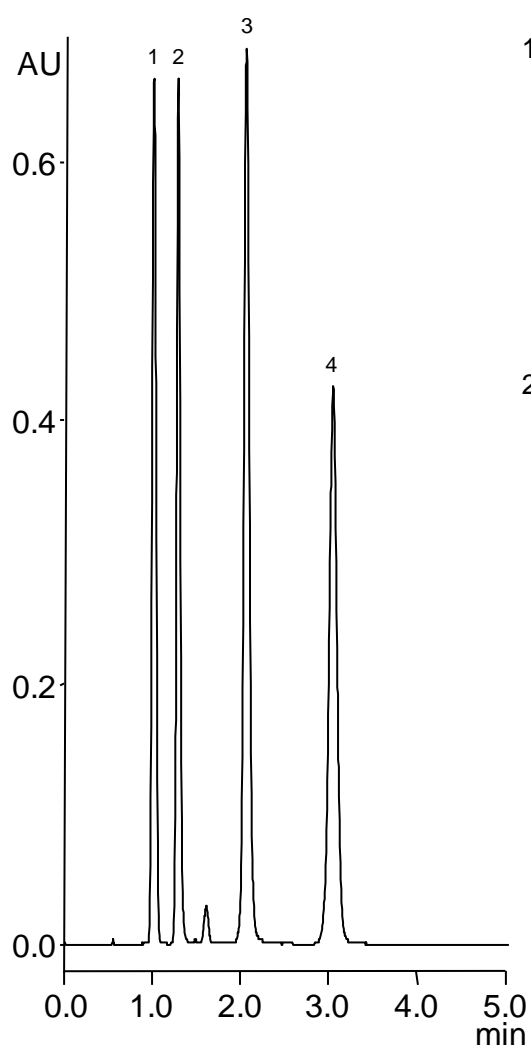


1.

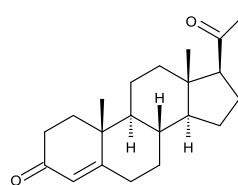


Progesterone

2.

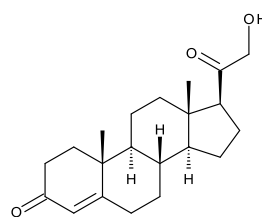


Deoxycorticosterone

3.

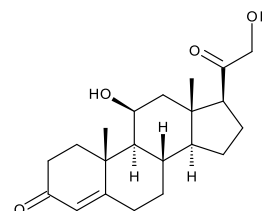

Corticosterone

4.

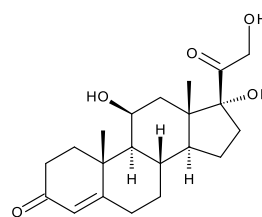

Hydrocortisone

Column : Alcyon SFC Triart Diol (5 μ m, 12 nm)
150 X 4.6 mm I.D.

Eluent : CO₂/ethanol (80/20)

Flow rate : 3.0 mL/min

Temperature : 40°C

Detection : UV at 254 nm

Back Pressure : 13.8 MPa (2000 psi)

Injection : 5 μ L (0.8 mg/mL)