



ステロイド
Steroids

X150803P

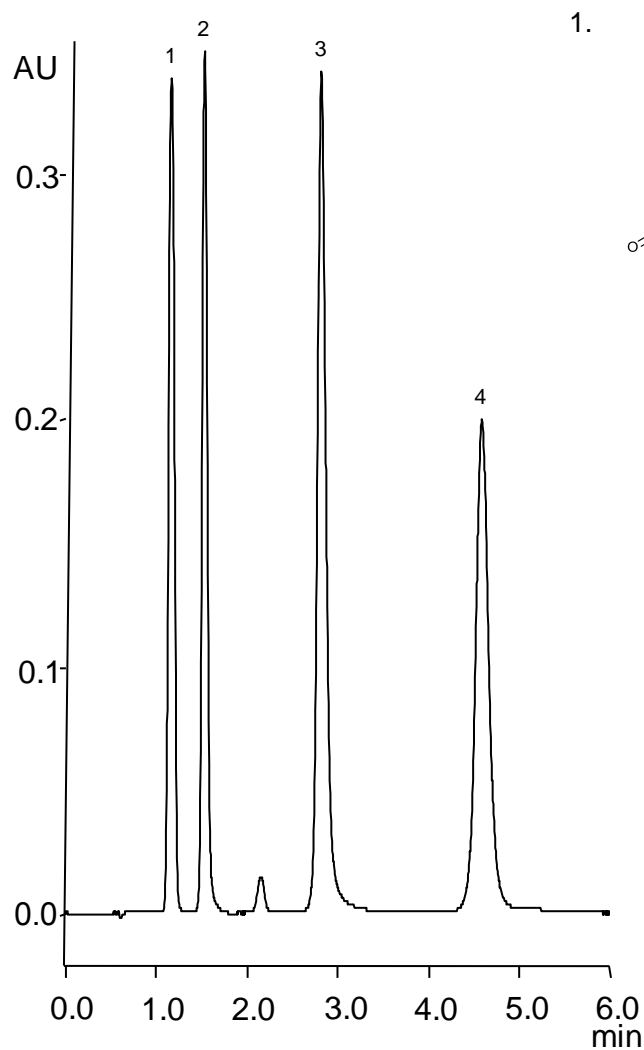


1.

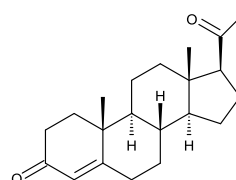


Progesterone

2.

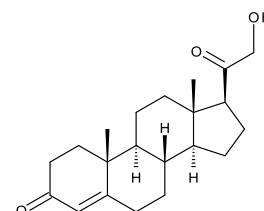


Deoxycorticosterone

3.

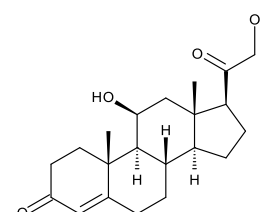

Corticosterone

4.

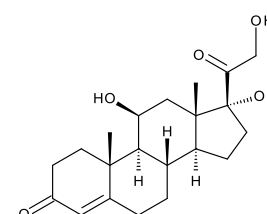

Hydrocortisone

| | |
|---------------|---|
| Column | : Alcyon SFC Triart Diol (5 μm, 12 nm) 150 X 4.6 mm I.D. |
| Eluent | : CO ₂ /2-propanol (80/20) |
| Flow rate | : 3.0 mL/min |
| Temperature | : 40°C |
| Detection | : UV at 254 nm |
| Back Pressure | : 13.8 MPa (2000 psi) |
| Injection | : 5 μL (0.8 mg/mL) |