


酸性非ステロイド性抗炎症薬

Acidic non-steroidal anti-inflammatory drugs

X151118F

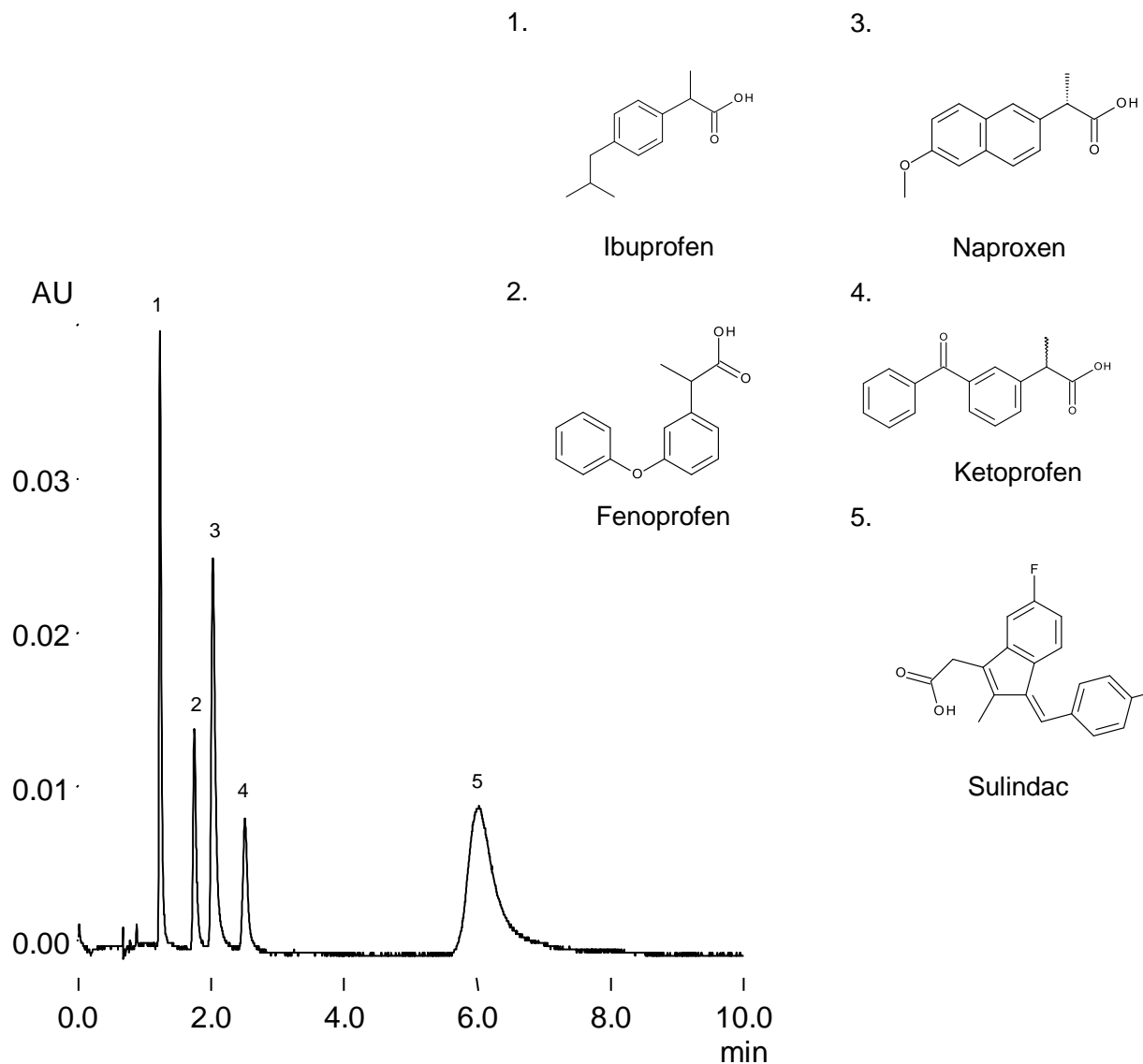

Column	: Alcyon SFC SIL (5 μm, 12 nm) 150 X 4.6 mm I.D.
Eluent	: CO <sub>2</sub> /methanol (0.1%DEA) (90/10)
Flow rate	: 3.0 mL/min
Temperature	: 40°C
Detection	: UV at 254 nm
Back Pressure	: 13.8 MPa (2000 psi)
Injection	: 1 μL (0.044mg/mL ~ 5.32 mg/mL)