


$\beta$ -遮断薬  
 $\beta$ -Blocker

X160420B

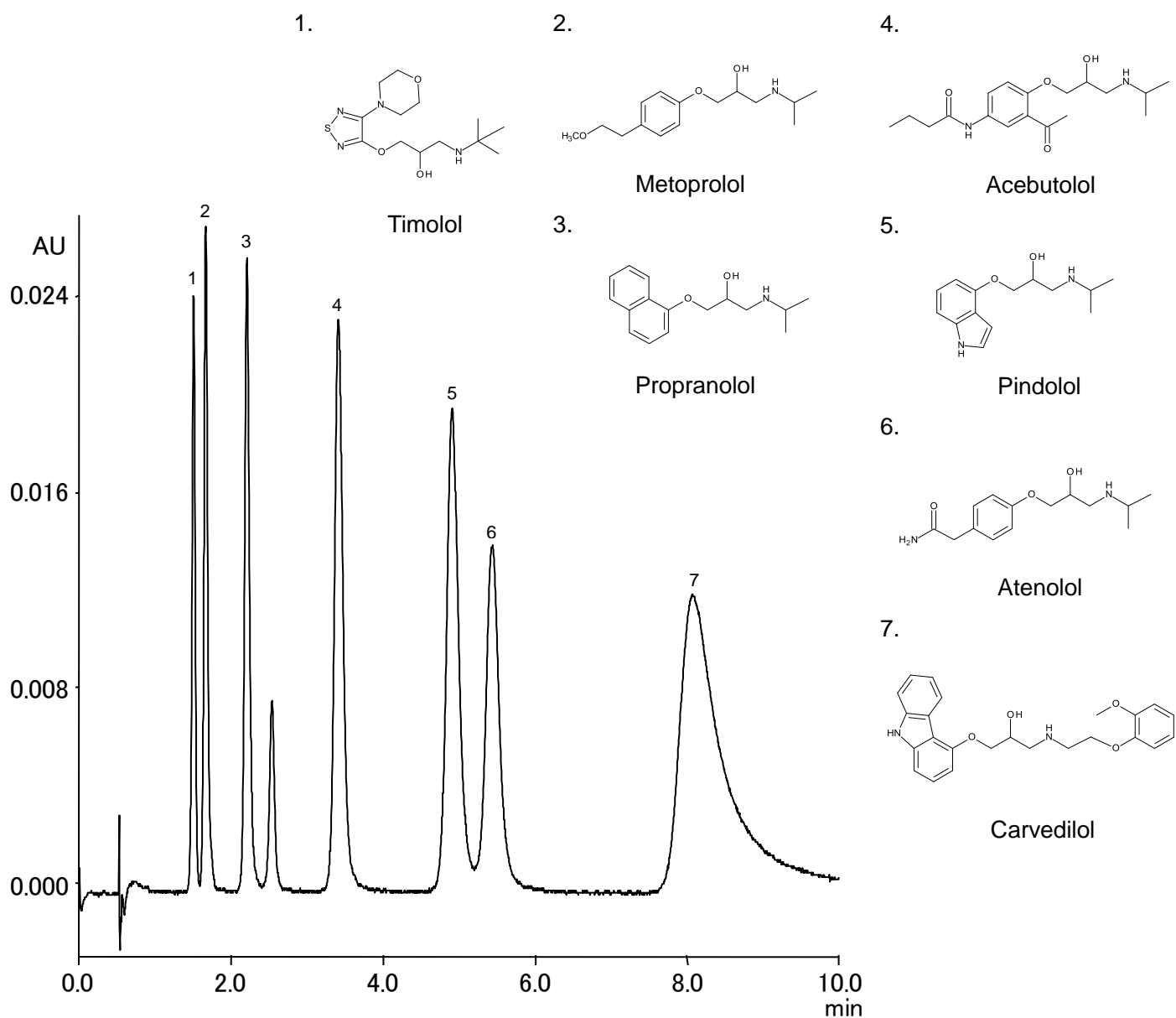

Column : Alcyon SFC Triart Diol (5  $\mu$ m, 12 nm)  
150 X 4.6 mm I.D.

Eluent : CO<sub>2</sub>/methanol 0.5%DEA (80/20)

Flow rate : 3.0 mL/min

Temperature : 40°C

Detection : UV at 254 nm

Back Pressure : 13.8 MPa (2000 psi)

Injection : 1  $\mu$ L (0.25 mg/mL ~ 4.6 mg/mL)