


DL-Phenyl ethyl amine

KR17010603

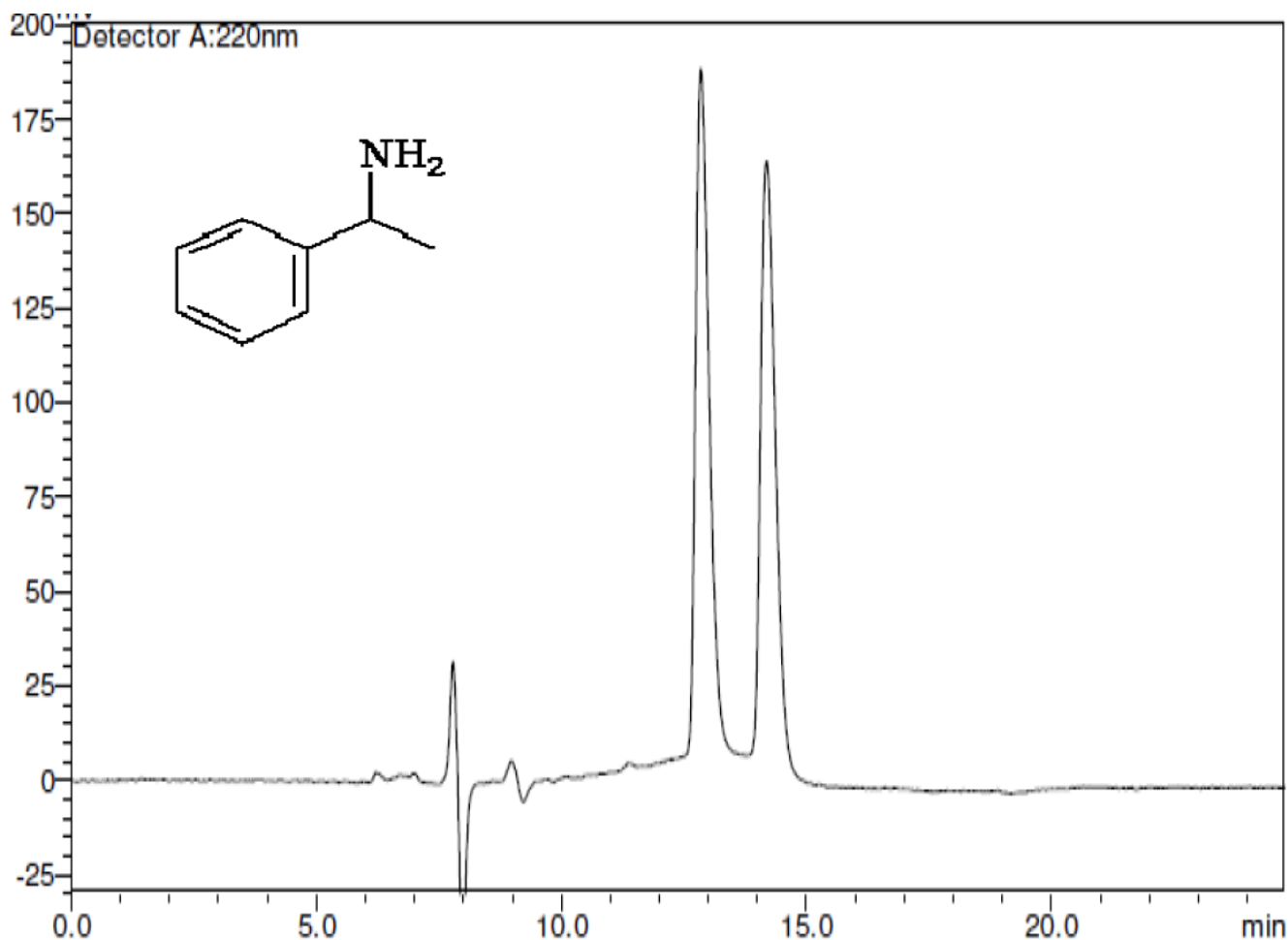

Retention Time (RT)	Resolution (Rs)	Name
12.845	-	Peak -1
14.183	2.50	Peak -2

Column : CHIRALART Cellulose-SB  
 250 X 4.6 mm I.D. 5 $\mu$ m  
 Eluent : n-hexane/Isopropyl alcohol/Diethylamine (90/10/0.2)  
 Flow rate : 0.5ml/min  
 Temperature : 25 $^{\circ}$ C  
 Detection : UV 220nm  
 Injection : 20 $\mu$ l (2.5mg/ml)