


Rosuvastatin

KR17010602

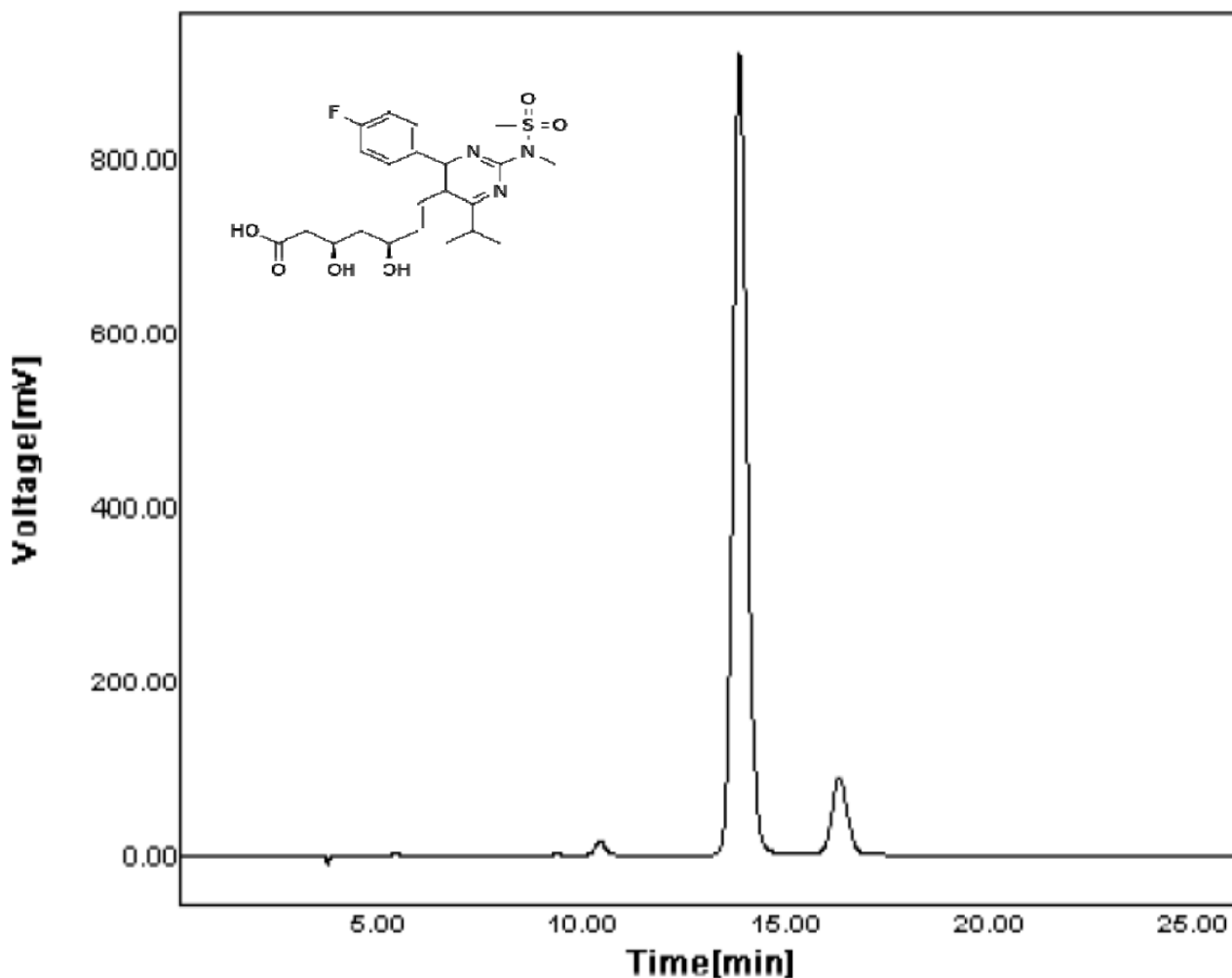

Retention Time (RT)	Resolution (Rs)	Name
10.383	-	-
13.800	6.52	Rosuvastatin
16.250	3.73	-

Column : CHIRALART Cellulose-SB  
 250 X 4.6 mm I.D. 5 $\mu$ m  
 Eluent : n-hexane/Ethanol/Trifluoroacetic acid (85/15/0.1)  
 Flow rate : 1 ml/min  
 Temperature : 25 $^{\circ}$ C  
 Detection : UV 242nm  
 Injection : 20 $\mu$ l (0.5mg/ml)