


Rivaroxaban

KR17010601

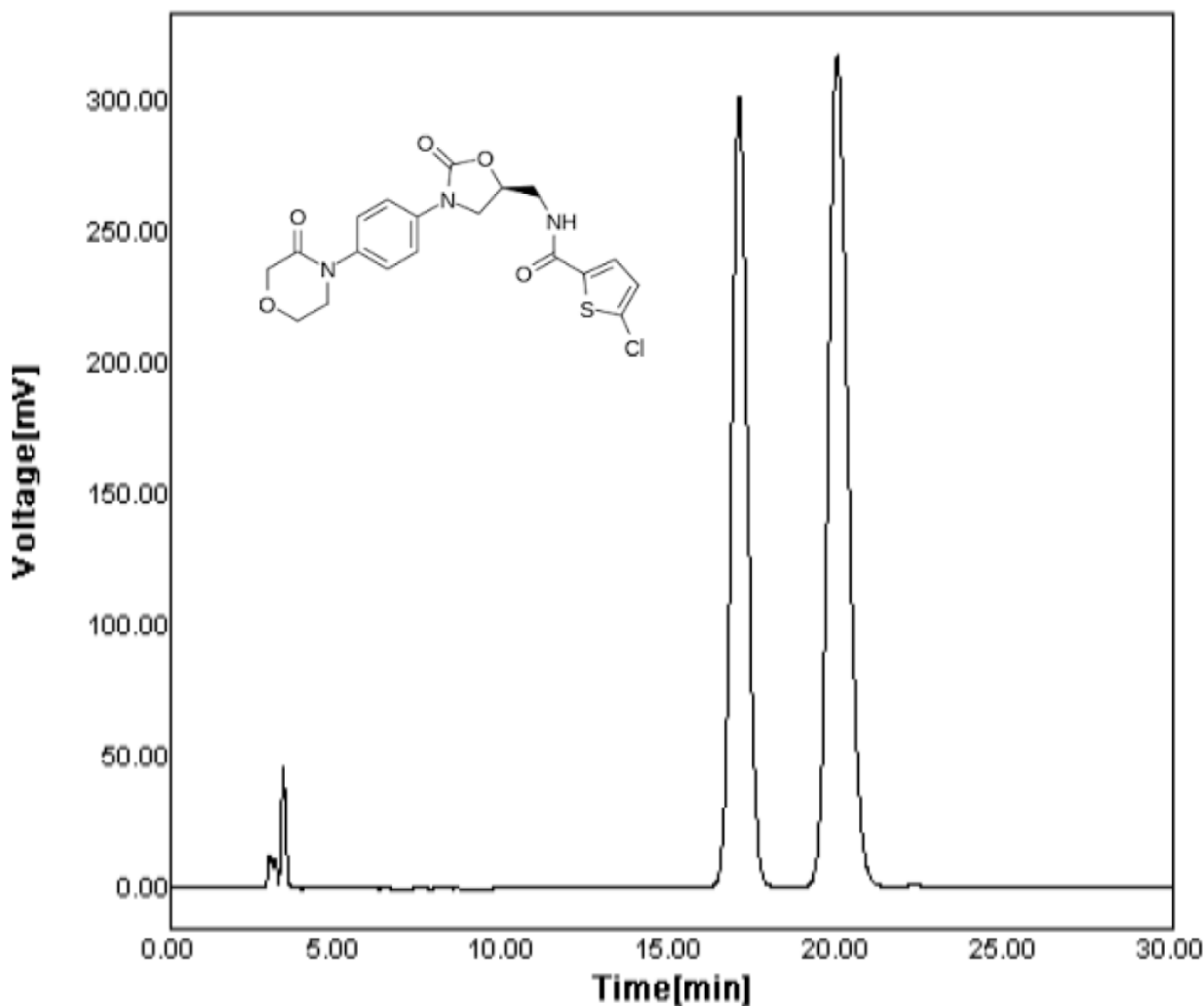

Retention Time (RT)	Resolution (Rs)	Name
17.050	-	RBN-III Isomer
19.967	2.99	RBN-III

Column : CHIRALART Cellulose-SB  
 250 X 4.6 mm I.D. 5 $\mu$ m  
 Eluent : n-hexane/Ethanol/Trifluoroacetic acid (50/50/0.1)  
 Flow rate : 1 ml/min  
 Temperature : 35 $^{\circ}$ C  
 Detection : UV 250nm  
 Injection : 20 $\mu$ l (0.5mg/ml)