


Edoxaban Tosylate Hydrate

KR17010607

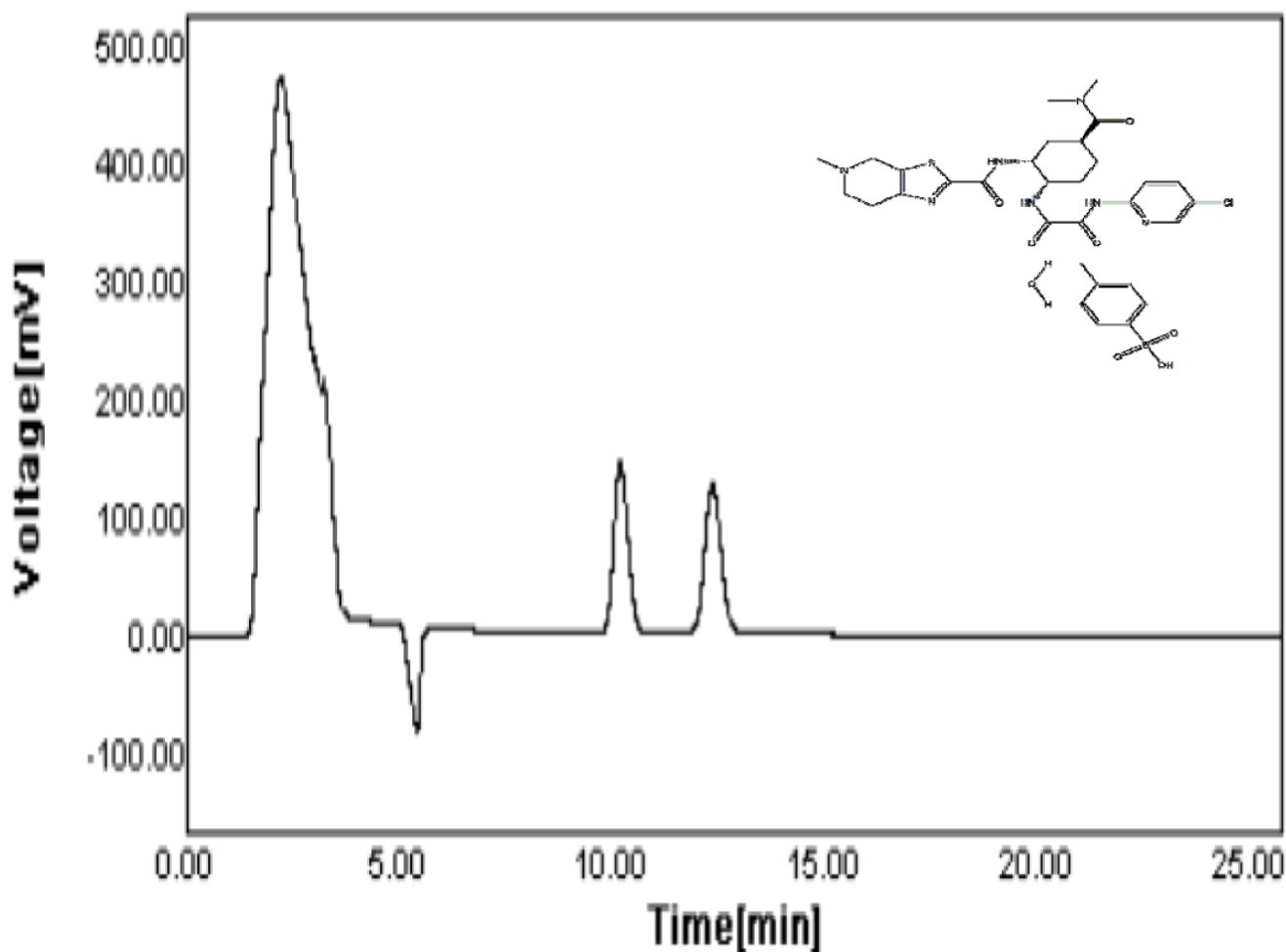

Retention Time (RT)	Resolution (R _s)	Name
10.187	-	Isomer
12.350	3.28	KSM-1

Column : CHIRALART Cellulose-SC
 250 X 4.6 mmI.D. 5 μ m
 Eluent : n-Hexane/Ethanol/Ethanolamine (55/45/0.01)
 Flow rate : 1.0ml/min
 Temperature : 25 $^{\circ}$ C
 Detection : UV 210nm
 Injection : 20 μ l (0.5mg/ml diluent :Methanol)