


Fluvastatin

KR17020101

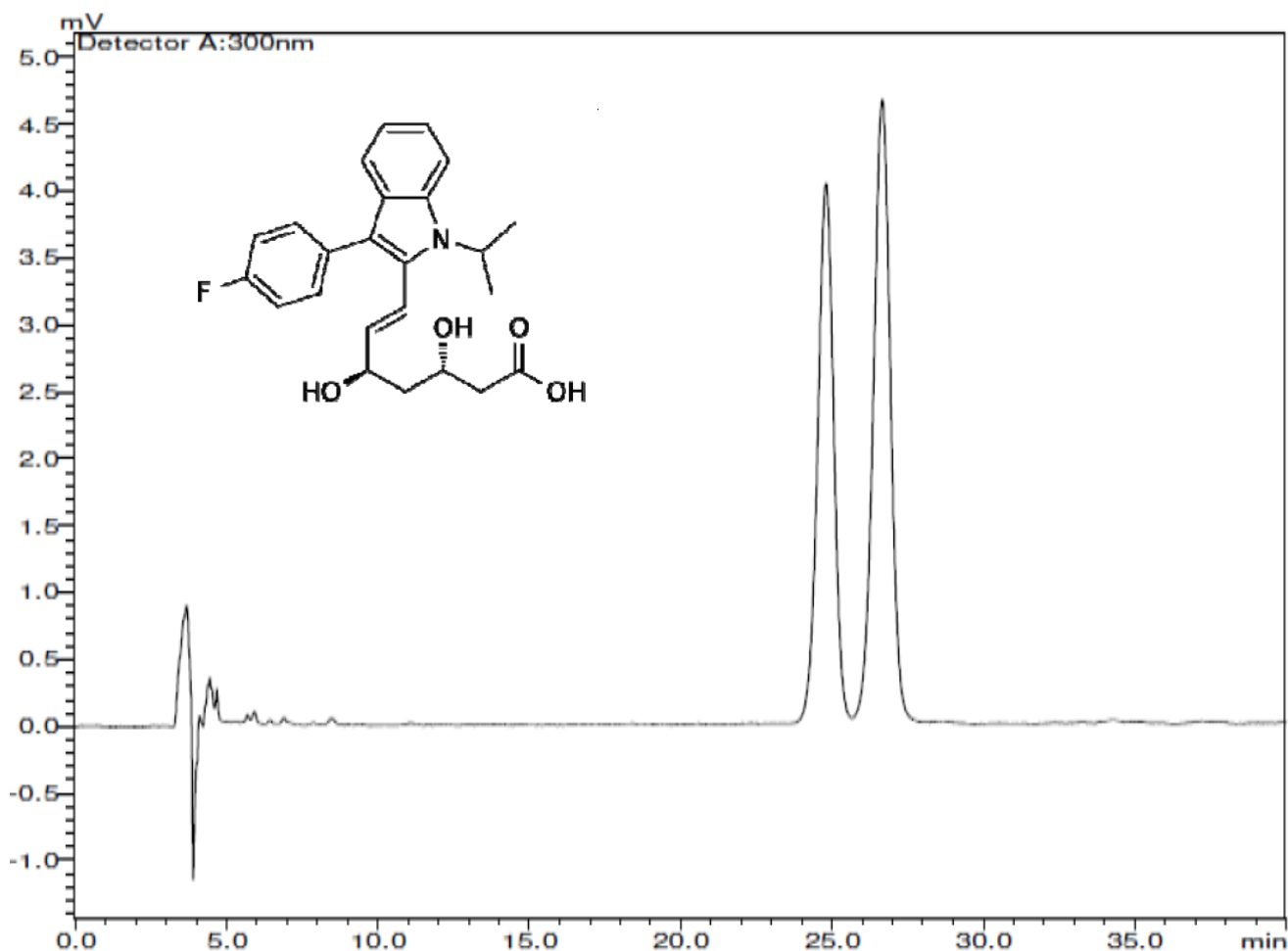

Retention Time (RT)	Resolution (Rs)	Name
24.798	-	3S 5R Fluvastatin
26.653	1.86	3R 5S Fluvastatin

Column : CHIRALART Cellulose-SB
250 X 4.6 mm I.D. 5 μ m

Eluent : 0.1% Formic acid in mixture of Water/Acetonitrile(65/35)

Flow rate : 1ml/min

Temperature : 25 $^{\circ}$ C

Detection : UV 300nm

Injection : 20 μ l (5mg/ml)