


Dorzolamide HCl

KR17020102

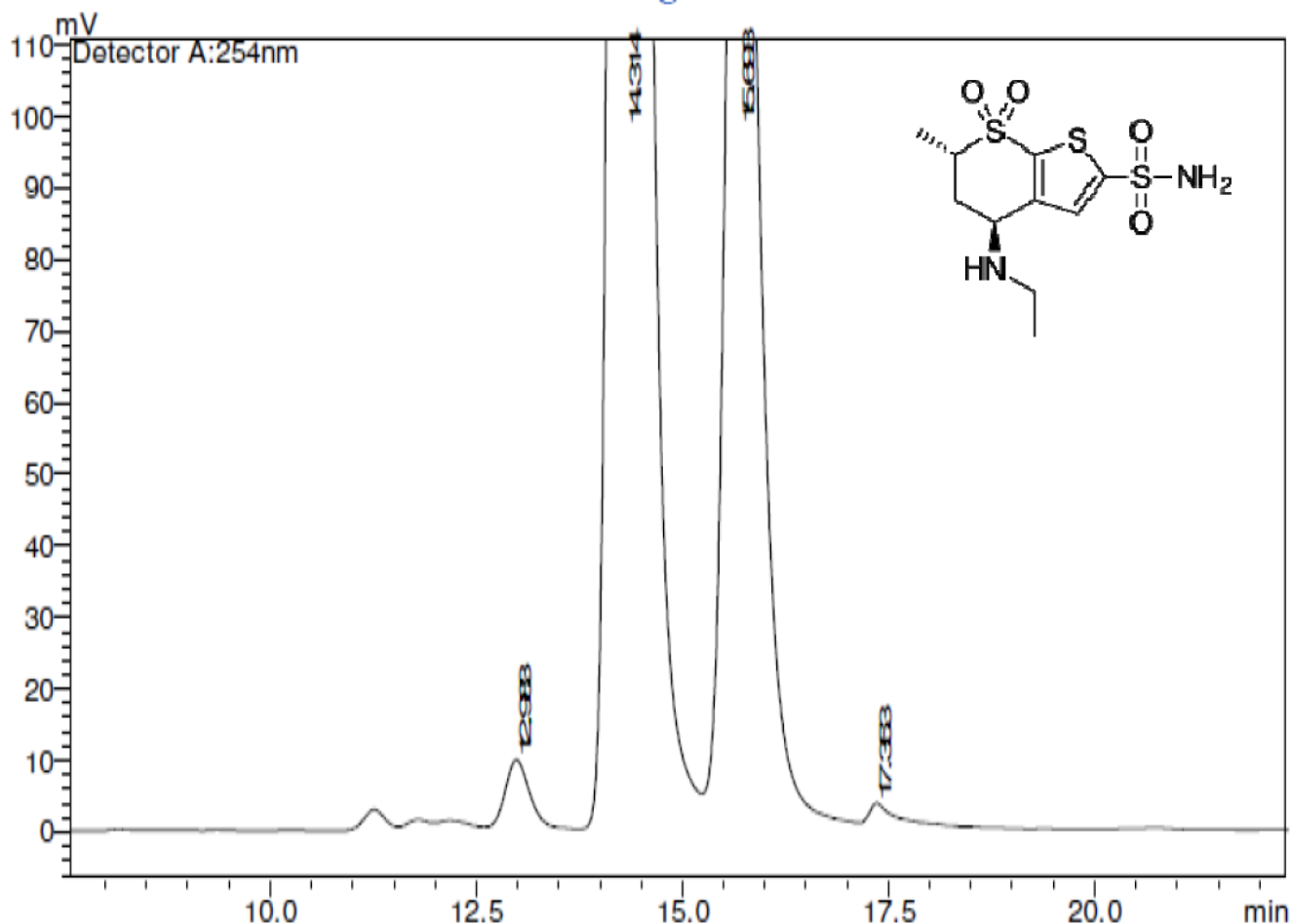

<i>Retention Time (RT)</i>	<i>Resolution (Rs)</i>	<i>Name</i>
12.983	-	Enantiomer-1
14.314	2.38	Dorzolamide
15.693	2.22	Enantiomer-3
17.353	3.10	Enantiomer-4

Column : CHIRALART Cellulose-SC
250 X 4.6 mm I.D. 5 μ m

Eluent : n-hexane/Ethanol/Diethylamine (80/20/0.1)

Flow rate : 1ml/min

Temperature : 25 $^{\circ}$ C

Detection : UV 254nm

Injection : 20 μ l (0.25mg/ml)