


ヌクレオチド  
Nucleotides


F161226

Metal free column

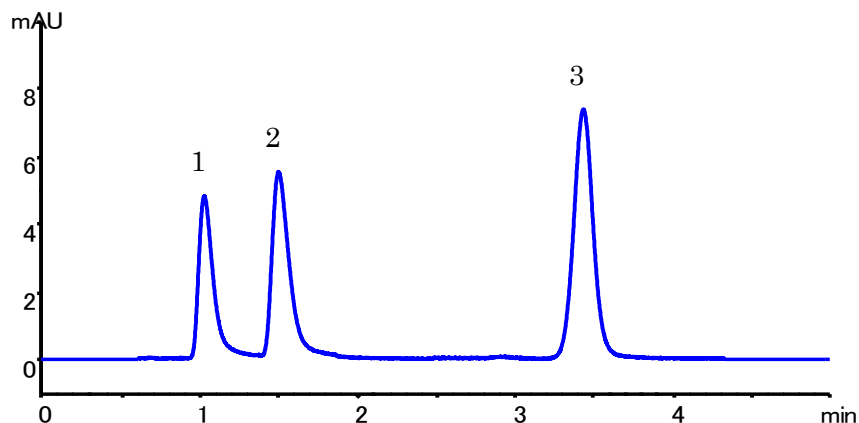

1. 5'-ATP
2. 5'-ADP
3. 5'-AMP

Standard column

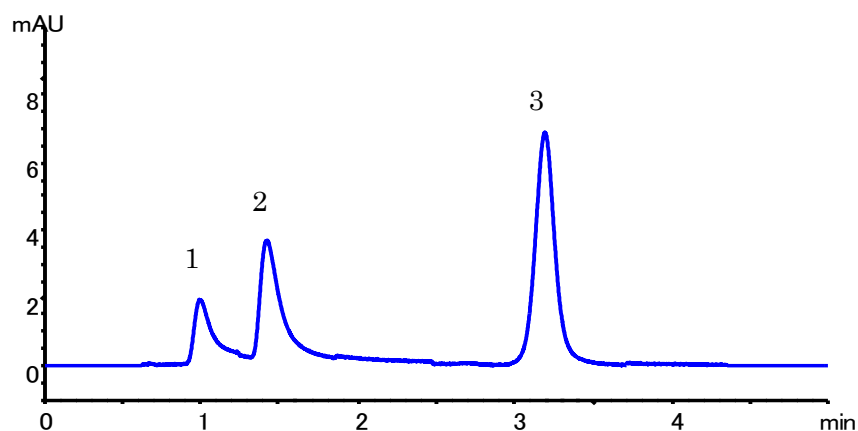

Column	: YMC-Triart C18 (3 $\mu$ m, 12 nm) 50 X 2.1 mm I.D.
Eluent	: 5 mM HCOONH <sub>4</sub>
Flow rate	: 0.21 mL/min
Temperature	: 25°C
Detection	: UV at 265 nm
Injection	: 1 $\mu$ L (10 $\mu$ g/mL)