


セチリジン (レボセチリジン)


Cetirizine (Levocetirizine)

F180214C

$\alpha$	1.1
Rs	2.4

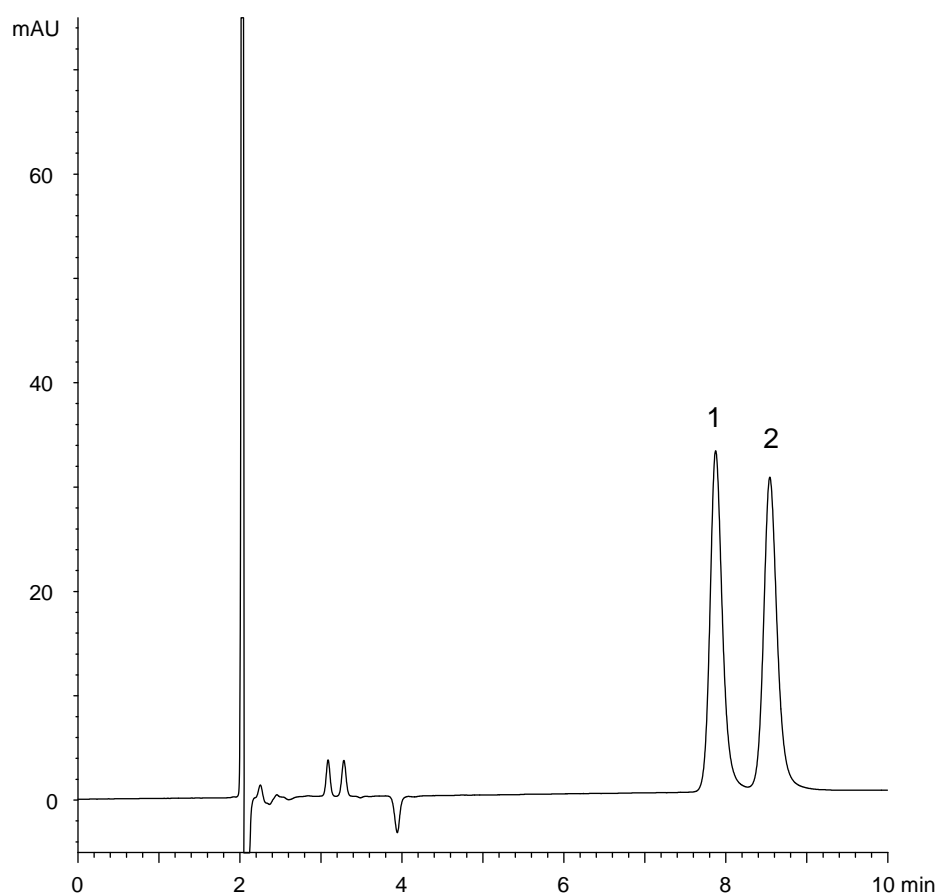


1

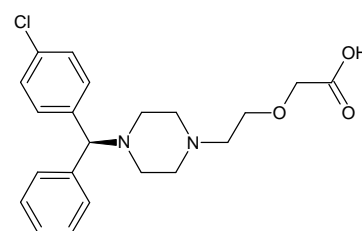

(S)-Cetirizine

2

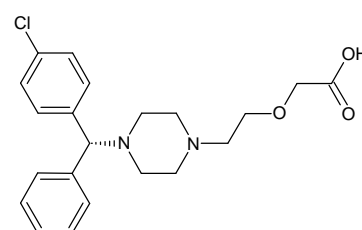

(R)-Cetirizine  
(Levocetirizine)

Column : CHIRAL ART Cellulose-SB (3  $\mu$ m)  
150 X 4.6 mmI.D.  
Eluent : acetonitrile/formic acid/diethylamine (100/0.1/0.1)  
Flow rate : 1.0 mL/min  
Temperature : 25°C  
Detection : UV at 230 nm  
Injection : 2  $\mu$ L (200  $\mu$ g/mL)